



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
 6-19-2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006053 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/19/14 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 5/22/14 and 6/19/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
 The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


6/19/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 19th, day of June, 2014.

Edma Escarcega

Notary Public

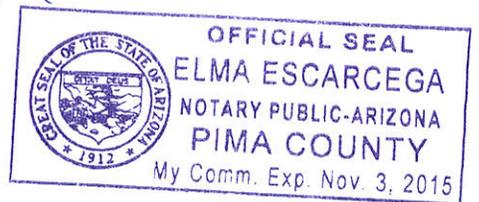

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
SC  
TS4  
6.19.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/19/14 TIME 829

INTOXILYZER SERIAL # 80-006053 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
06/19/2014

08:31:59

2 of 7 SC  
TS4

---

Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:32:49
Diagnostic Test	Pass	08:33:24
Air Blank	0.000	08:33:51
0.100 Cal Check	0.097	08:34:12
Air Blank	0.000	08:34:41
Subject Test	0.000	08:35:06
Air Blank	0.000	08:35:33
Wait	ABT*	
Air Blank	0.000	08:36:29

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
06/19/2014

08:36:53

3 of 7 sc  
TS4

---

Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:37:45
Diagnostic Test	Pass	08:38:20
Air Blank	0.000	08:38:47
0.100 Cal Check	0.098	08:39:08
Air Blank	0.000	08:39:38
Subject Test	RFI*	08:40:07
Air Blank	IPS**	08:40:19

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
06/19/2014

08:44:28

4 of 7  
sc  
154

---

Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	08:45:23
Diagnostic Test	Pass	08:45:58
Air Blank	0.000	08:46:25
0.100 Cal Check	0.098	08:46:46
Air Blank	0.000	08:47:16
Subject Test	0.108*	08:50:23
Air Blank	0.000	08:51:02
Air Blank	0.000	08:51:29
Subject Test	XXX**	08:52:02
Air Blank	0.000	08:52:35

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
06/19/2014 08:53:00

5 of 7  
sc  
TS4

-----  
Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:53:29
0.100 Cal Check	0.096	08:53:50
Air Blank	0.000	08:54:19

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7<sub>sc</sub>  
TS4  
6-19-2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/19/14 TIME 829  
Intoxilyzer Serial # 80-006053 Location DUT HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
06/19/2014

08:29:10

7 of 7 sc  
By

-----  
Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:29:39
0.100 Cal Check	0.097	08:30:00
Air Blank	0.000	08:30:29

Preventative Maintenance  
Performed -- Timer Reset