



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TSy  
9.5.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006053 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/14/14 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 7/17/14 and 8/14/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


8/19/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 19th, day of August, 20 14.

Elma Escarcega

Notary Public

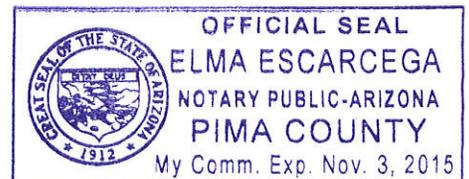

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7 sc  
TS4  
9.5.2014

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 8/14/14 TIME 1045

INTOXILYZER SERIAL # 80-006053 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
08/14/2014

10:49:44

2 of 7 sc  
BH

---

Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:50:34
Diagnostic Test	Pass	10:51:09
Air Blank	0.000	10:51:36
0.100 Cal Check	0.097	10:51:57
Air Blank	0.000	10:52:27
Subject Test	0.000	10:53:25
Air Blank	0.000	10:53:52
Wait	ABT*	
Air Blank	0.000	10:54:26

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
08/14/2014

10:58:41

3 of 7 sc  
BH

---

Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:59:46
Diagnostic Test	Pass	11:00:21
Air Blank	0.000	11:00:48
0.100 Cal Check	0.096	11:01:09
Air Blank	0.000	11:01:39
Subject Test	RFI*	11:01:55
Air Blank	IPS**	11:02:02

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
08/14/2014

11:03:13

4 of 7<sub>sc</sub>  
Tsh

---

Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:03:55
Diagnostic Test	Pass	11:04:30
Air Blank	0.000	11:04:57
0.100 Cal Check	0.098	11:05:18
Air Blank	0.000	11:05:48
Subject Test	0.133*	11:08:55
Air Blank	0.000	11:09:27
Air Blank	0.000	11:09:54
Subject Test	XXX**	11:10:18
Air Blank	0.000	11:10:52

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
08/14/2014 11:12:48

5 of 7  
sk  
TS4

-----  
Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:13:16
0.100 Cal Check	0.095	11:13:37
Air Blank	0.000	11:14:06

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 SC  
TSY  
9.5.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 8/14/14 TIME 1045  
Intoxilyzer Serial # 80-006053 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006053  
Core Version: 8105.48  
08/14/2014 10:45:17

7 of 7 SC  
By

-----  
Standard Lot#: 25913100A2  
Last Changed By: A. MARQUIS #47674  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:45:46
0.100 Cal Check	0.096	10:46:07
Air Blank	0.000	10:46:36

Preventative Maintenance  
Performed -- Timer Reset