



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006053 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

WLL  
6/29/16

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-7-2016 Cal 1 N/A AC Cal 2 N/A AC

Between the dates of \_\_\_\_\_ and \_\_\_\_\_ the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 6-7-2016.  
The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


6-29-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 29<sup>th</sup>, day of June, 2016.

Elma Escarcega

Notary Public

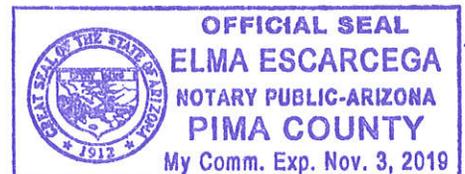

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 3  
TSH

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

WLL  
6/29/16

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28844 AGENCY Tucson Police Department  
DATE 6-7-2016 TIME N/A  
Intoxilyzer Serial # 80-006053 Location CrimeLab

- (  ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (  ) 2. Intoxilyzer 8000 display reads "~~PUSH BUTTON TO START~~" *"Fault #18 error, Contact CMI, Inc"*
- ( ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ) 4. Air blank completed.
- ( ) 5. Calibration check completed. Test results 0. \_\_\_\_\_ AC.
- ( ) 6. Air blank completed.
- ( ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

2 of 3  
TS4

INTOXILYZER SERVICE RECORD

WLL  
6/29/16

Date: 6-7-2016

QAS: Terry Gallegos 28847

Intoxilyzer Serial Number: 80-006053

Reason For Service: "Fault #18 error, Contact CMI, Inc."  
on display when turned on. (EEPROM fail, DSP fail,  
Analytical fail).

Date Out of Service: 6-7-2016

Calibration Check completed prior to removal from service

Date: NO

SQAP completed prior to removal from service

Date: NO

Corrective Action:  Sent to CMI, Inc.  
Records downloaded

**CMI** INC.  
**Service Evaluation Form**

3 of 3  
T84  
WLL 6/29/16

This form MUST be completed and enclosed with instrument to be serviced.

Note: Please ship items in their original shipping container.

**Contact information:** Customer Number 105075 (contact Customer Service)  
Name TERRY GALLEGOS Phone: (520) 837-7677  
Fax: (520) 791-4642 Email: TERRY.GALLEGOS@TUCSONAZ.GOV

**Bill to Address:** FINANCE/ACCOUNTS PAYABLE  
PO BOX 27450  
TUCSON, AZ 85726-7450

**Ship to Address:** TUCSON POLICE DEPARTMENT CRIMELAB  
1306 W. MIRAGE MILE  
TUCSON, AZ 85705

Instrument Serial Number: 80-006053


Detailed Description of Problem:  
"Fault #18 error, Contact CMI, Inc." on display when  
turned on.

\*\*\*Hazardous Material Warning! - DO NOT return gas cylinders with instrument! \*\*\*

Note: An evaluation fee of \$39.50 for handheld instruments or \$79.00 for IR evidential instruments will apply to estimates which are not repaired.

I Authorize Repairs Up To:  All  \$250  \$500  \$750  Other \$ \_\_\_\_\_  
Purchase Order Number 37768

**Authorized By:**

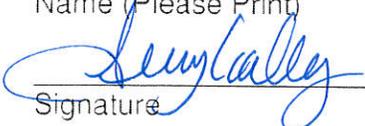TERRY GALLEGOS TOXICOLOGY COORDINATOR Ship to:  
Name (Please Print) Title  
 6.22.2016  
Signature Date

CMI, Inc.  
Attn: Service Dept.  
316 East Ninth Street  
Owensboro, KY 42303

No, please send estimate before repairs are made. An estimate will be faxed before performing any repairs and may cause delays in service.