



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
8.4.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006054 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/17/14 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 6/19/14 and 7/17/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


7/18/14

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 18th, day of July, 2014.

Elma Escarcega

Notary Public

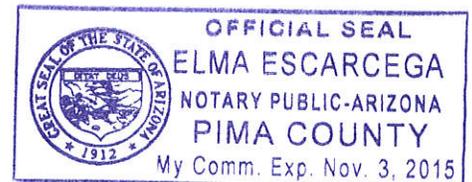

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
sc  
TSY  
8.4.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 7/17/14 TIME 1043  
INTOXILYZER SERIAL # 80-006054 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/17/2014

10:44:35

2 of 7<sub>sc</sub>  
186

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:45:22
Diagnostic Test	Pass	10:45:57
Air Blank	0.000	10:46:25
0.100 Cal Check	0.098	10:46:45
Air Blank	0.000	10:47:15
Subject Test	0.000	10:47:29
Air Blank	0.000	10:47:56
Wait	ABT*	
Air Blank	0.000	10:50:09

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/17/2014

10:50:36

3 of 7  
sc  
TS4

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:51:45
Diagnostic Test	Pass	10:52:20
Air Blank	0.000	10:52:48
0.100 Cal Check	0.098	10:53:08
Air Blank	0.000	10:53:38
Subject Test	RFI*	10:54:18
Air Blank	IPS**	10:54:31

\*RFI Detect

\*\*Improper Sample

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/17/2014

10:55:29

4 of 7  
SC  
B6

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:56:19
Diagnostic Test	Pass	10:56:54
Air Blank	0.000	10:57:21
0.100 Cal Check	0.097	10:57:42
Air Blank	0.000	10:58:12
Subject Test	0.135*	11:01:19
Air Blank	0.000	11:01:51
Air Blank	0.000	11:02:19
Subject Test	XXX**	11:02:46
Air Blank	0.000	11:03:20

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/17/2014 11:03:43

5 of 7  
SC  
B4

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:04:12
0.100 Cal Check	0.097	11:04:33
Air Blank	0.000	11:05:02

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077  
SC  
BS4  
8.4.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 7/17/14 TIME 1043  
Intoxilyzer Serial # 80-006054 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/17/2014 10:43:09

7 of 7 sc  
TD4

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:43:37
0.100 Cal Check	0.098	10:43:58
Air Blank	0.000	10:44:28

Preventative Maintenance  
Performed -- Timer Reset