



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

*Old*  
7/30/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-006054 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7.16.2015 Cal 1 0.098 AC Cal 2 0.099 AC

Between the dates of 6.18.2015 and 7.16.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

*Terry Gallegos*

Terry S. Gallegos


7.30.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 30<sup>th</sup>, day of July, 2015.

*Elma Escarcega*

Notary Public

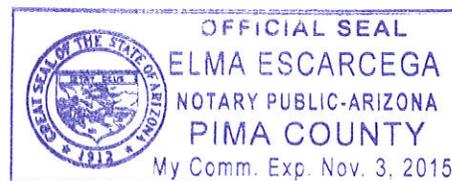

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 9  
TSh  
7/20/15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 7-16-2015 TIME 1221

INTOXILYZER SERIAL # 80-006054 LOCATION DUI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---




---



---



---

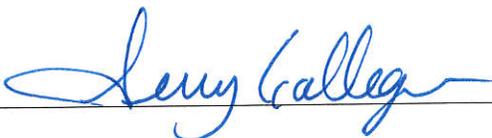SIGNATURE 

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/16/2015

12:23:45

2 of 9  
TSh  
dls

---

Standard Lot#: 25414100A4  
Last Changed By: C. MORIN #49542

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:24:44
Diagnostic Test	Pass	12:25:19
Air Blank	0.000	12:25:46
0.100 Cal Check	0.099	12:26:07
Air Blank	0.000	12:26:36
Subject Test	0.000	12:26:51
Air Blank	0.000	12:27:18
Wait	ABT*	
Air Blank	0.000	12:27:52

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/16/2015

12:28:33

30f9  
TB4  
JLB

---

Standard Lot#: 25414100A4  
Last Changed By: C. MORIN #49542

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:29:26
Diagnostic Test	Pass	12:30:01
Air Blank	0.000	12:30:28
0.100 Cal Check	RFI*	12:30:46
Air Blank	IPS**	12:30:59

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/16/2015 12:31:49

4 of 9  
754  
DLS

Standard Lot#: 25414100A4  
Last Changed By: C. MORIN #49542

QAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:32:38
Diagnostic Test	Pass	12:33:13
Air Blank	0.000	12:33:40
0.100 Cal Check	0.099	12:34:01
Air Blank	0.000	12:34:31
Subject Test	0.069*	12:37:38
Air Blank	0.000	12:38:10
Air Blank	0.000	12:38:37
Subject Test	XXX**	12:38:53
Air Blank	0.000	12:39:27

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/16/2015 12:40:30

50f9  
B4  
DD

-----  
Standard Lot#: 25414100A4  
Last Changed By: C. MORIN #49542  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:40:58
0.100 Cal Check	0.099	12:41:19
Air Blank	0.000	12:41:49
Air Blank	0.000	12:44:17

\*Input Timed Out

*maintenance updated B4*

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 9  
TS4  
AD  
7/30/15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 7-16-2015 TIME 1221  
Intoxilyzer Serial # 80-006054 Location DUI-HQ

- (  ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (  ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (  ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (  ) 4. Air blank completed.
- (  ) 5. Calibration check completed. Test results 0. 098 AC.
- (  ) 6. Air blank completed.
- (  ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/16/2015 12:21:45

*Fof 9  
TS4  
[Signature]*

-----  
Standard Lot#: 25414100A4  
Last Changed By: C. MORIN #49542  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:22:17
0.100 Cal Check	0.098	12:22:37
Air Blank	0.000	12:23:07

Preventative Maintenance  
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8 of 9  
754  
AB  
7/30/15

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 07 04 15  
 QAS: MORIN 49542  
 Intoxilyzer Serial Number: 80-006054  
 New Lot Number: 25414100A4  
 New Expiration Date: 10/5/2016

- () 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- () 2. Remove the old gas tank and install the new one.
- () 3. Hit the "ESC" key twice and swipe the QAS card.
- () 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
 Select "D" for dry gas and enter the following information:  
     Enter Target Value: 0.100  
     Enter Control Lot #: From the tank label  
     Enter Exp Date: From the tank label (mm/dd/yyyy)  
     Edit info Y/N? Y or N  
     Update Std Info? Y ENTER  
     "Please wait saving..."
- () 5. Hit the "ESC" key twice.
- () 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
 Select "D" for dry control test. After selection is made press ENTER.
- () 7. Air blank completed.
- () 8. Calibration check completed. Test results 0. 099 AC.
- () 9. Air blank completed.
- () 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

[Signature] 49542

90f9  
TS4  
Add

INTOXILYZER 8000  
Instrument Initialization  
20:24:11 07/04/2015

OLD TRAVE  
C. MORIN 49542

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/04/2015 20:55:16

Standard Lot#: ~~25414100A4~~ 11 514100A2  
Last Changed By: C. MORIN #49542

QAS: C. MORIN #49542  
TUCSON PD

Operator: C. MORIN #49542  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	20:55:49
0.100 Cal Check	0.099	20:56:10
Air Blank	0.000	20:56:39

Preventative Maintenance  
Performed -- Timer Reset

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-006054  
Core Version: 8105.48  
07/04/2015 21:00:16

Standard Lot#: 25414100A4  
Last Changed By: C. MORIN #49542

QAS: C. MORIN #49542  
TUCSON PD

Operator: C. MORIN #49542  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	21:00:47
0.100 Cal Check	0.099	21:01:08
Air Blank	0.000	21:01:37

C. MORIN 49542

Preventative Maintenance  
Performed -- Timer Reset