



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-006054 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

WLL
6/29/16

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-17-2016 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 5-20-2016 and 6-17-2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


6-29-2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 29th, day of June, 2016.

Elma Escarcega

Notary Public

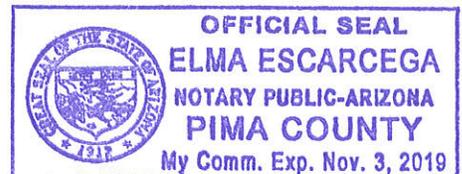

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10710
TS4
WLL
6/29/16

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6-17-2016 TIME 1155
INTOXILYZER SERIAL # 80-006054 LOCATION Crime Lab

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006054
Core Version: 8105.48
06/17/2016

11:58:25

2 of 10
B4
WLL

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:59:12
Diagnostic Test	Pass	11:59:47
Air Blank	0.000	12:00:15
0.100 Cal Check	0.102	12:00:35
Air Blank	0.000	12:01:05
Subject Test	0.000	12:01:32
Air Blank	0.000	12:02:00
Wait	ABT*	
Air Blank	0.000	12:02:42

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006054
Core Version: 8105.48
06/17/2016

12:03:14

3 of 10
TS4
WLL

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:04:12
Diagnostic Test	Pass	12:04:47
Air Blank	0.000	12:05:14
0.100 Cal Check	RFI*	12:05:30
Air Blank	IPS**	12:05:38

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006054
Core Version: 8105.48
06/17/2016

12:08:42

4 of 10
TS6
WLL

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:09:28
Diagnostic Test	Pass	12:10:03
Air Blank	0.000	12:10:30
0.100 Cal Check	0.101	12:10:51
Air Blank	0.000	12:11:21
Subject Test	0.169*	12:14:28
Air Blank	0.000	12:15:04
Air Blank	0.000	12:15:32
Subject Test	XXX**	12:16:01
Air Blank	0.000	12:16:35

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006054
Core Version: 8105.48
06/17/2016 12:16:45

Sofio
TS
WLL

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:17:14
0.100 Cal Check	0.100	12:17:35
Air Blank	0.000	12:18:04

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

60F10
TS4
WLL
6/29/16

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6-17-2016 TIME 1155
Intoxilyzer Serial # 80-006054 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-006054
Core Version: 8105.48
06/17/2016 11:55:17

*7 of 10
TS4
WLL*

Standard Lot#: 25414100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:55:48
0.100 Cal Check	0.101	11:56:09
Air Blank	0.000	11:56:39

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

80f10
B4

WLL
6/29/16

Date: 6-7-2016

QAS: Terry Gallegos 28847

Intoxilyzer Serial Number: 80-006054

Reason For Service: Notified by note that the instrument was not working. I checked the calibration (G-2) and a subject test performed. Both were good tests. Officer Estes had three tests in which voltage/current failed Diagnostics. The battery in the vehicle will be checked for proper voltage.

Date Out of Service: Not out of service.

Calibration Check completed prior to removal from service

Date: _____

SQAP completed prior to removal from service

Date: _____

Corrective Action: Sent to CMI, Inc.

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

9 of 10
TS4
WLL
6/29/16

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6-7-2016 TIME 1250
Intoxilyzer Serial # 80-006054 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000
 Instrument Initialization
 12:23:10 06/06/2016

INTOXILYZER 8000
 Location: TUCSON PD
 Serial Number: 80-006054
 Core Version: 8105.48
 06/07/2016 12:50:56

Standard Lot#: 25414100A4
 Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
 TUCSON PD

Operator: T. GALLEGOS #28847
 TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:51:26
0.100 Cal Check	0.100	12:51:46
Air Blank	0.000	12:52:16

Preventative Maintenance NOT Performed -- Timer NOT Reset

INTOXILYZER 8000
 Location: TUCSON PD
 Serial Number: 80-006054
 Core Version: 8105.48
 06/07/2016 12:52:37

Standard Lot#: 25414100A4
 Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
 TUCSON PD

Operator: T. GALLEGOS #28847
 TUCSON PD

Subject: TEST,SQA,TEST
 DOB: 11/11/1987
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:53:37
Diagnostic Test	Pass	12:54:14
Air Blank	0.000	12:54:42
0.100 Cal Check	0.100	12:55:02
Air Blank	0.000	12:55:32
Subject Test	0.000	12:55:58
Air Blank	0.000	12:56:25
Wait	ABT*	
Air Blank	0.000	12:57:01

*Sequence Aborted

Not a Successfully
 Completed
 Test Sequence

10 of 10
 TSh
 WLL