



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

[Handwritten signature]

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 64-001226 (#8) maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 9-8-2006 Cal 1 0.080 AC Cal 2 0.079 AC

Date: 10-6-2006 Cal 1 0.082 AC Cal 2 0.080 AC

Date: — Cal 1 — AC Cal 2 — AC

Between the dates of 9-8-2006 and 10-6-2006 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Terry Gallegos

Terry S. Gallegos

10-30-2006

ADPS QAS Permit # QAS 976

SUBSCRIBED AND SWORN to before me this 30 day of October, 2006.

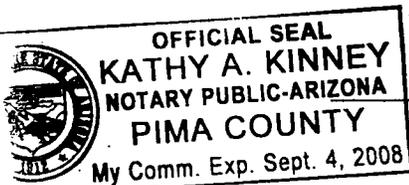

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10/6
Dy

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD QUALITY ASSURANCE PROCEDURE *ju*

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 10.6.2006 TIME 1550
INTOXILYZER SERIAL # 64-001226 (#8) LOCATION DUI HQ

(✓) 1. Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- (✓) 1. DVM Test check. Setting should be between .010 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .355.
- (✓) 2. Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- (✓) 3. Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- * (✓) 4. Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- (✓) 5. Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject Test result 0.000 AC
- (✓) 2. Error recognition logic system functioning
Invalid test printed
- (✓) 3. Proper sample recognition system
Invalid sample printed ✓
Deficient sample printed ✓
- (✓) 4. Standard Calibration Check standard 0.080 AC
Results: 0.082 AC and 0.080 AC

Instrument operating accurately and properly. YES NO

COMMENTS: * 1616 → 1621

SIGNATURE *Terry Gallegos*

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DPS Form Exh E-4 (Rev. 05-1)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

Tucson Police Department Crime Laboratory

3076
Bj

FEDERAL SIGNAL CORP. / CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 64-001226
10/06/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	15:55
CAL. CHECK	.082	15:55
AIR BLANK	.000	15:55

#8

SUBJECT'S NAME

DUI #9

TIME FIRST OBSERVED

INSTRUMENT LOCATION

Gallegos 20047

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER™ INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

FEDERAL SIGNAL CORP. / CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 64-001226
10/06/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	15:56
CAL. CHECK	.080	15:57
AIR BLANK	.000	15:57

#8

SUBJECT'S NAME

DUI #9

TIME FIRST OBSERVED

INSTRUMENT LOCATION

Gallegos 20047

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER™ INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

Hoff
D4
file

SN 64-001226
E47.02
INVALID TEST
SAMPLE INTRODUCED
AT IMPROPER TIME.

10/06/06
16:15

#8

SUBJECT'S NAME

DUI HQ
INSTRUMENT LOCATION

TIME FIRST OBSERVED

Gallegos 23047
OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER[®] INSTRUMENT PRINTER CARD



© 1986 by CMI INC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

5066
BH

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE

file

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 10-6-2006 TIME 1550
INTOXILYZER SERIAL # 64-001226 (#8) LOCATION DUI HQ

- (✓) 1. Pour a standard alcohol concentration solution into a clean dry simulator jar and assemble the simulator. Ensure that a tight seal has been made. Standard value: 0. 080 AC.
- (✓) 2. Turn on the simulator and allow the temperature to reach 34°C ± 0.2°C.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓) 6. Push Start Test button or press enter on keyboard.
- (✓) 7. Insert card in response to display.
- (✓) 8. Air blank completed.
- (✓) 9. Standard Calibration Check completed. Test results 0. 082 AC and 0. 080 AC.
- (✓) 10. Air blank completed.
- (✓) 11. When test complete, remove printed record. Attach the record to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE Terry Gallegos

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DPS Form Exh E-2 (Rev. 05-1)

FEDERAL SIGNAL CORP. / CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 64-001226
10/06/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	16:12
*SUBJECT TEST	.075	16:14
AIR BLANK	.000	16:14

* DEFICIENT SAMPLE - VALUE
PRINTED WHO HIGHEST OBTAINED.

#8

SUBJECT'S NAME

DUI HQ

INSTRUMENT LOCATION

TIME FIRST OBSERVED

Gallegos 28847

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER™ INSTRUMENT PRINTER CARD

© 1986 by CMI INC

FEDERAL SIGNAL CORP. / CMI INC
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 64-001226
10/06/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	16:10
INVALID SAMPLE	.XXX	16:11
AIR BLANK	.000	16:11

#8

SUBJECT'S NAME

DUI HQ

INSTRUMENT LOCATION

TIME FIRST OBSERVED

Gallegos 28847

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet breath



INTOXILYZER™ INSTRUMENT PRINTER CARD

© 1986 by CMI INC

6076 D4
jfe