



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

ph

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 5000 Serial Number 66-001759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 2/23/06 Cal 1 0.079 AC Cal 2 0.078 AC

Date: 3/23/06 Cal 1 0.084 AC Cal 2 0.084 AC

Date: N/A Cal 1                      AC Cal 2                      AC

Between the dates of 2/23/06 and 3/23/06 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on 3/23/06.  
The instrument was returned to service on N/A.

Jennifer L. Coombs  
Jennifer L. Coombs

3/30/06

ADPS QAS Permit # QAS 29886

SUBSCRIBED AND SWORN to before me this 30, day of March, 2006.

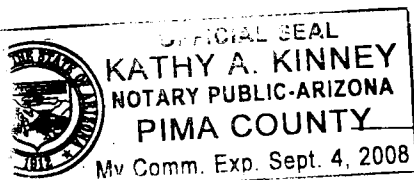

Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R9-14-404.A.6

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

## EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 5000\*

134

AGENCY Tucson Police QA SPECIALIST Jennifer Coombs 761160  
(Print Name)INTOXILYZER SERIAL # 66-CC1759 LOCATION DUI HQDATE March 23 19 2006 TIME 1615

☒ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

## DIAGNOSTIC TESTS

- ☒ DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .356.
- ☒ Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- ☒ Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- ☒ Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- ☒ Date check. Mode selection switch S11 on or keyboard menu selection "E".

## OPERATIONAL TESTS

- ☒ Alcohol free subject Test result 0.000 AC
- ☒ Error recognition logic system functioning  
Invalid test printed
- ☒ Proper sample recognition system  
Invalid sample printed ✓  
Deficient sample printed ✓
- N/A Completeness of sample purge with collection tube\*\*
- ☒ Calibration standard 0.680 AC
- Results: 0.084 AC & 0.084 AC
- Instrument operating properly and accurately. YES ☒ NO ☐

COMMENTS: \* Moved clock 1612 → 1615SIGNATURE Jennifer L. Coombs

0. N/A alcohol concentration standard collected for subsequent analysis.\*\*

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD  
\*\*THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR  
SAMPLE CAPTURE.

DHS/DSLS/Form C117 (Rev. 7-93)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

TEST	BREQ VALUE	TIME
TEST BLANK	.000	16:22
SUBJECT TEST	.000	16:22
AIR BLANK	.000	16:23

$$DM = .3536$$

SUBJECT'S NAME \_\_\_\_\_

Dull HQ

TIME FIRST OBSERVED

46160  
OPERATOR

[illegible]

**EMI** INC.

'INTOXILYZER' INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

[illegible]

SUBJECT'S NAME:

Dull HD

TIME FIRST OBSERVED

46160 OPERATOR

[illegible]

INTOXILYZER<sup>®</sup> INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

# Tucson Police Department Crime Laboratory *TSy*

*3 of 6 file*

TUCSON POLICE DEPT.  
INTOXILYZER - ALCOHOL ANALYZER /  
AZ MODEL 5000 SN 66-001759  
03/23/06 16:21 /

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:20
CAL. CHECK	.084	16:21
AIR BLANK	.000	16:21

SUBJECT'S NAME  
*DUI HQ*

INSTRUMENT LOCATION

TIME FIRST OBSERVED

*461160*  
OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC

TUCSON POLICE DEPT.  
INTOXILYZER - ALCOHOL ANALYZER /  
AZ MODEL 5000 SN 66-001759  
03/23/06 16:20 /

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:19
CAL. CHECK	.084	16:19
AIR BLANK	.000	16:20

SUBJECT'S NAME

*DUI HQ*  
INSTRUMENT LOCATION

TIME FIRST OBSERVED

*461160*  
OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC

4 of 6  
file  
104

SN 66-001759  
E767.14  
INVALID TEST

03/23/86/  
16:24/

SAMPLE INTRODUCED  
AT IMPROPER TIME.

SUBJECT'S NAME	
TIME FIRST OBSERVED	INSTRUMENT LOCATION DUI HQ
OPERATOR 461160	
ADDITIONAL INFORMATION AND/OR REMARKS	

INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC

**CMI** INC

EXHIBIT P

5 of 6  
file

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 5000\*

BSH

STANDARD CALIBRATION CHECK PROCEDURE

Agency Tucson Police Department Date 3/23/06 Time 11:15  
Intoxilyzer Serial # 666-001759 Location DUI HQ  
QA Specialist Jennifer Combs 461160  
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach 34°C ± .2°C.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- ✓ ✓ 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- ✓ ✓ 6. Push Start Test button or press enter on keyboard.
- ✓ ✓ 7. Insert card in response to display.
- ✓ ✓ 8. Air blank completed.
- ✓ ✓ 9. Calibration check completed. Test results 0.084AC & 0.084 AC
- ✓ ✓ 10. Air blank completed.
- ✓ ✓ 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Jennifer L. Combs

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

TUCSON POLICE DEPT.  
INTOXILYZER - ALCOHOL ANALYZER /  
AZ MODEL 5000 SN 66-001759  
03/23/06 16:26 /

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:24
INVALID SAMPLE	.000	16:26
AIR BLANK	.000	16:26

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME

DUI HQ

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet Breath



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC

TUCSON POLICE DEPT.  
INTOXILYZER - ALCOHOL ANALYZER /  
AZ MODEL 5000 SN 66-001759  
03/23/06 16:31 /

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:28
*SUBJECT TEST	.105	16:30
AIR BLANK	.000	16:31

\* DEFICIENT SAMPLE - VALUE  
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME

DUI HQ

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

134

Date: March 23, 2006

QAS: Jennifer Coombs 461160

Intoxilyzer Serial Number: 66-001759

Reason For Service: Calibration checks have increased.

Date Out of Service: 3/23/06

Calibration Check completed prior to removal from service

Date: 3/23/06

SQAP completed prior to removal from service

Date: 3/23/06

Corrective Action:

☒ Sent to Applied Electronics