



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

fh

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 5000 Serial Number 66-004014 (#17) maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 3/23/06 Cal 1 0.080 AC Cal 2 0.079 AC

Date: 4/20/06 Cal 1 0.080 AC Cal 2 0.078 AC

Date: N/A Cal 1                      AC Cal 2                      AC

Between the dates of 3/23/06 and 4/20/06 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Jennifer L. Coombs  
Jennifer L. Coombs

4/21/06

ADPS QAS Permit # QAS 29886

SUBSCRIBED AND SWORN to before me this 21, day of April, 2006.

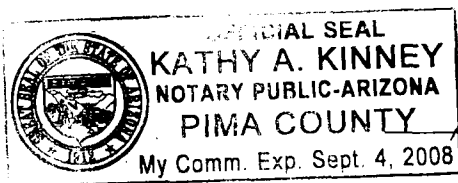

Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R9-14-404.A.6

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

## ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 5000\*AGENCY Tucson Police QA SPECIALIST Jennifer Coombs 46160  
(Print Name)INTOXILYZER SERIAL # 66-004014 (#17) LOCATION DUI HQ  
DATE April 20 19 2006 TIME 1449

☒ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

## DIAGNOSTIC TESTS

☒ DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .507.

☒ Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".

☒ Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".

☒ Clock time check. Mode selection switch S10 on or keyboard menu selection "E".

☒ Date check. Mode selection switch S11 on or keyboard menu selection "E".

## OPERATIONAL TESTS

☒ Alcohol free subject Test result 0.000 AC

☒ Error recognition logic system functioning  
Invalid test printed

☒ Proper sample recognition system  
Invalid sample printed  
Deficient sample printed

N/A Completeness of sample purge with collection tube\*\*

☒ Calibration standard 0.180 AC

Results: 0.050 AC & 0.078 AC

Instrument operating properly and accurately. YES ☒ NO ☐

COMMENTS: \* moved clock 1446 → 1449

SIGNATURE

Jennifer L Coombs  
0. N/A alcohol concentration standard collected for subsequent analysis.\*\*

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD  
\*\*THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR  
SAMPLE CAPTURE.

DHS/DSLS/Form C117 (Rev. 7-93)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

# Tucson Police Department Crime Laboratory

3 of 6  
file

#17

CMI INC  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-004014  
04/20/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	15:00
CAL. CHECK	.078	15:01
AIR BLANK	.000	15:01

SUBJECT'S NAME

DUI HO

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

CMI

INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

#17

CMI INC  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-004014  
04/20/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	14:59
CAL. CHECK	.080	14:59
AIR BLANK	.000	15:00

SUBJECT'S NAME

DUI HO

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

CMI

INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

4 of 6  
JH

#17

SN 66-004014  
E647.05  
INVALID TEST  
SAMPLE INTRODUCED  
AT IMPROPER TIME.

04/20/06  
14:55

SUBJECT'S NAME	
TIME FIRST OBSERVED	DUI HQ INSTRUMENT LOCATION
46160 OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	

INTOXILYZER<sup>®</sup> INSTRUMENT PRINTER CARD



© 1986 by CMI INC

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 5000\*

STANDARD CALIBRATION CHECK PROCEDURE

Agency Tucson Police Department Date 4/20/06 Time 1449  
Intoxilyzer Serial # 66-00404 (#17) Location DUI #76  
QA Specialist Jennifer Coombs 46160  
(Print Name)

- ( ✓ ) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- ( ✓ ) 2. Turn on the simulator and allow the temperature to reach  $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$ .
- ( ✓ ) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- ( ✓ ) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- ✓ ✓ 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- ✓ ✓ 6. Push Start Test button or press enter on keyboard.
- ✓ ✓ 7. Insert card in response to display.
- ✓ ✓ 8. Air blank completed.
- ✓ ✓ 9. Calibration check completed. Test results 0.080 AC & 0.078 AC
- ✓ ✓ 10. Air blank completed.
- ✓ ✓ 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- ( ✓ ) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Jennifer L. Coombs

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

#17  
CMI INC  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-004014  
04/20/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	15:03
INVALID SAMPLE	.XXX	15:04
AIR BLANK	.000	15:04

TIME FIRST OBSERVED  
SUBJECT'S NAME  
INSTRUMENT LOCATION

40160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet Breath

INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC

CMI INC

#17  
CMI INC  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-004014  
04/20/06

TEST	BAC VALUE	TIME
AIR BLANK	.000	15:05
*SUBJECT TEST	.073	15:07
AIR BLANK	.000	15:07

\* DEFICIENT SAMPLE - VALUE  
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME

TIME FIRST OBSERVED

INSTRUMENT LOCATION

40160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC

CMI INC

6 of 6  
gle