



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

BM

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-004015 (#18) maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 12/20/05 Cal 1 0.077 AC Cal 2 0.078 AC

Date: 1/4/06 Cal 1 0.079 AC Cal 2 0.078 AC

Date: N/A Cal 1 AC Cal 2 AC

Between the dates of 12/20/05 and 1/4/06
the instrument was operating properly and accurately.


The instrument was taken out of service for repair on 1/4/06.
The instrument was returned to service on N/A.

Jennifer L. Coombs
Jennifer L. Coombs

1/12/06

ADPS QAS Permit # QAS 29886

SUBSCRIBED AND SWORN to before me this 12, day of
January, 2006.

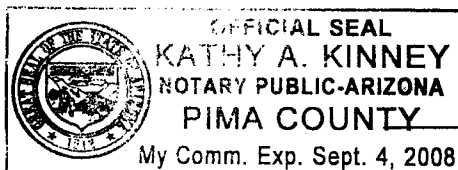

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R9-14-404.A.6

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES *SH*
INTOXILYZER MODEL 5000*AGENCY Tucson Police QA SPECIALIST Jennifer Coombs 46160
(Print Name)INTOXILYZER SERIAL # 66-004015 (#18) LOCATION LABDATE January 4 18 fe 2006 TIME 1530☒ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

☒ DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .319.☒ Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".☒ Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".☒ Clock time check. Mode selection switch S10 on or keyboard menu selection "E".☒ Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

☒ Alcohol free subject Test result 0.000 AC☒ Error recognition logic system functioning
Invalid test printed☒ Proper sample recognition system
Invalid sample printed
Deficient sample printed

N/A Completeness of sample purge with collection tube**

☒ Calibration standard 0.080 ACResults: 0.079 AC & 0.078 ACInstrument operating properly and accurately. YES ☒ NO ☐COMMENTS: * rotated year dial to 2006
** printer not printing clearlySIGNATURE Jennifer L. Coombs

0. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD
**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR
SAMPLE CAPTURE.

DHS/DSLS/Form C117 (Rev. 7-93)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 3
fle

2 of 3
file

INSTRUMENT #18	
DATE	
TIME	
LOCATION	
OPERATOR	
SUBJECT'S NAME	
TIME FIRST OBSERVED	
INSTRUMENT LOCATION	
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	
SN 666-004015	
46160	
LTB	
INTOXILYZER® INSTRUMENT PRINTER CARD	
© 1986 by CMI INC.	

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE

Agency Tucson Police Department Date 1/4/06 Time 1530
Intoxilyzer Serial # 66-004015(*18) Location LAB
QA Specialist Jennifer Coombs 46160
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- ✓ ✓ 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- ✓ ✓ 6. Push Start Test button or press enter on keyboard.
- ✓ ✓ 7. Insert card in response to display.
- ✓ ✓ 8. Air blank completed.
- ✓ ✓ 9. Calibration check completed. Test results 0.079 AC & 0.078 AC
- ✓ ✓ 10. Air blank completed.
- ✓ ✓ * 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Jennifer L. Coombs

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

* Printer not printing clearly



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

TSY

Date: January 3, 2006

QAS: Jennifer Coombs 46660

Intoxilyzer Serial Number: 66-004015 (#18)

Reason For Service: On 12/31/05 Sgt Skeenes 34967 emailed Terry Gallegos that there was a printer problem with Intoxilyzer 5000 #18 at ODW. He was instructed to put an out of service sign on it. I picked up the instrument on 1/3/06 and brought it to the lab.

Date Out of Service: 1/4/06

Calibration Check completed prior to removal from service

Date: 1/4/06

SQAP completed prior to removal from service

Date: 1/4/06

Corrective Action:

☒ Sent to Applied Electronics
