



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-004798 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program. 124

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 7/7/05 Cal 1 0.079 AC Cal 2 0.078 AC

Date: 8/8/05 Cal 1 0.077 AC Cal 2 0.077 AC

Date: — Cal 1 — AC Cal 2 — AC

Between the dates of 7/7/05 and 8/8/05 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.


Irene Doran

Irene Doran

8/9/05

ADPS QAS Permit # QAS 29586

SUBSCRIBED AND SWORN to before me this 9, day of August, 2005.

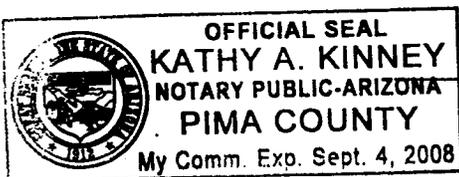

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R9-14-404.A.6

EXHIBIT Q

1086
2x

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES *154*
INTOXILYZER MODEL 5000*

AGENCY Tucson Police QA SPECIALIST Irene Doran 46331
(Print Name)

INTOXILYZER SERIAL # 66-004798 LOCATION ODD

DATE August 8th 2005 TIME 0935

Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .294.
- Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- Alcohol free subject Test result 0.000 AC
 - Error recognition logic system functioning
Invalid test printed
 - Proper sample recognition system
Invalid sample printed ✓
Deficient sample printed ✓
 - N/A Completeness of sample purge with collection tube**
 - Calibration standard 0.080 AC
Results: 0.077AC + 0.077AC
- Instrument operating properly and accurately. YES NO

COMMENTS: * 32 days since last SQAP.

SIGNATURE Irene Doran

0. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD
**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR SAMPLE CAPTURE.

DHS/DSLS/Form C117 (Rev. 7-93)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3076
PK

Tucson Police Department Crime Laboratory

TUCSON POLICE DEPARTMENT
 INTOXILYZER - ALCOHOL ANALYZER
 AZ MODEL 5000 SN 66-004798
 08/08/05 09:50

TEST BRAC VALUE TIME
 AIR BLANK .000 09:50
 CAL. CHECK .077 09:50
 AIR BLANK .000 09:50

SUBJECT'S NAME _____

TIME FIRST OBSERVED _____ INSTRUMENT LOCATION _____

OPERATOR **46331**

ADDITIONAL INFORMATION AND/OR REMARKS _____

INTOXILYZER® INSTRUMENT PRINTER CARD
 © 1986 by CMI INC.



TUCSON POLICE DEPARTMENT
 INTOXILYZER - ALCOHOL ANALYZER
 AZ MODEL 5000 SN 66-004798
 08/08/05 09:49

TEST BRAC VALUE TIME
 AIR BLANK .000 09:48
 CAL. CHECK .077 09:49
 AIR BLANK .000 09:49

SUBJECT'S NAME _____

TIME FIRST OBSERVED _____ INSTRUMENT LOCATION _____

OPERATOR **46331**

ADDITIONAL INFORMATION AND/OR REMARKS _____

INTOXILYZER® INSTRUMENT PRINTER CARD
 © 1986 by CMI INC.



496
PK

SN 66-004798
E767.14
INVALID TEST

08/08/05
09:41

SAMPLE INTRODUCED
AT IMPROPER TIME.

SUBJECT'S NAME

TIME FIRST OBSERVED

INSTRUMENT LOCATION

44331
OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI

596
EK

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

BM

STANDARD CALIBRATION CHECK PROCEDURE

Agency Tucson Police Department Date 8/8/05 Time 0935
Intoxilyzer Serial # 66-004798 Location ODD
QA Specialist Irene Doran 46331
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach 34°C+.2°C.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓) 6. Push Start Test button or press enter on keyboard.
- (✓) 7. Insert card in response to display.
- (✓) 8. Air blank completed.
- (✓) 9. Calibration check completed. Test results 0.077 + 0.077 AC
- (✓) 10. Air blank completed.
- (✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE Irene Doran

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

696
LA

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004798
08/08/05 09:46

TEST	BRAC VALUE	TIME
AIR BLANK	.000	09:43
*SUBJECT TEST	.130	09:45
AIR BLANK	.000	09:46

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME _____

TIME FIRST OBSERVED _____

OPERATOR 46331

INSTRUMENT LOCATION ODD

ADDITIONAL INFORMATION AND/OR REMARKS _____

INTOXILYZER® INSTRUMENT PRINTER CARD
© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004798
08/08/05 09:42

TEST	BRAC VALUE	TIME
AIR BLANK	.000	09:41
INVALID SAMPLE	.MOA	09:42
AIR BLANK	.000	09:42

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME _____

TIME FIRST OBSERVED _____

OPERATOR 46331

INSTRUMENT LOCATION ODD

ADDITIONAL INFORMATION AND/OR REMARKS Sweet breath

INTOXILYZER® INSTRUMENT PRINTER CARD
© 1986 by CMI INC.