



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

AM

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-004798 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 10/26/05 Cal 1 0.081 AC Cal 2 0.081 AC

Date: 11/23/05 Cal 1 0.078 AC Cal 2 0.079 AC

Date: — Cal 1 — AC Cal 2 — AC

Between the dates of 10/26/05 and 11/23/05 the instrument was working properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Irene Doran Johnson


Irene Doran Johnson

12/1/05

ADPS QAS Permit # QAS 29586

SUBSCRIBED AND SWORN to before me this 1, day of
December, 2005.

Kathy A Kinney
Notary Public

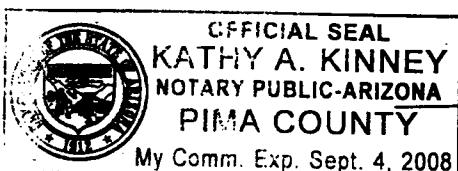

*Includes documents specified in A.A.C. R9-14-404.A.6

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

196
IJ

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES *AM*
INTOXILYZER MODEL 5000*

AGENCY Tucson Police QA SPECIALIST Irene D. Johnson 46331
(Print Name)

INTOXILYZER SERIAL # 66-004798 LOCATION ODD
DATE November 23 18 2005 TIME 0730

IJ

Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .317.
- Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- Alcohol free subject Test result 0.000 AC
- Error recognition logic system functioning
 - Invalid test printed
- Proper sample recognition system
 - Invalid sample printed
 - Deficient sample printed
- N/A Completeness of sample purge with collection tube**
 - Calibration standard 0.090 AC
 - Results: 0.078AC + 0.079AC
- Instrument operating properly and accurately. YES NO

COMMENTS:

SIGNATURE

Irene D. Johnson 46331

0. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR
SAMPLE CAPTURE.

TUCSON POLICE DEPARTMENT
 INTOXILYZER ALCOHOL ANALYZER
 AZ MODEL 5900 SN 66-04798
 11/23/05 07:48

TEST	BRAC VALUE	TIME
AIR BLANK	.000	07:47
SUBJECT TEST	.000	07:48
AIR BLANK	.000	07:48

SN 66-04798 11/23/05 07:49
 E767.14 EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE
 EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE
 EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE
 EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE
 EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE
 EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE
 EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE
 EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE EEEEEE
 Q307 -----
 Q306 -----
 Q305 -----
 Q302 -----
 Q303 -----
 Q304 -----
 Q301 -----

DVM = .317

29
HJ

SUBJECT'S NAME	ODD
INSTRUMENT LOCATION	TIME FIRST OBSERVED
OPERATOR	ADDITIONAL INFORMATION AND/OR REMARKS
46331	

SUBJECT'S NAME	ODD
INSTRUMENT LOCATION	TIME FIRST OBSERVED
OPERATOR	ADDITIONAL INFORMATION AND/OR REMARKS
46331	

CMI

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

396
IJ
E

Tucson Police Department Crime Laboratory

TUCSON POLICE DEPARTMENT INTOXILYZER - ALCOHOL ANALYZER		
AZ MODEL 5000	SN 66-004798	11/23/05
TEST	BRAC VALUE	TIME
AIR BLANK	.000	07:54
CAL. CHECK	.078	07:54
AIR BLANK	.000	07:55

TEST	BRAC VALUE	TIME
AIR BLANK	.000	07:54
CAL. CHECK	.078	07:54
AIR BLANK	.000	07:55

CMI

INTOXILYZER™ INSTRUMENT PRINTER CARD
© 1986 by CMI INC.

CMI

INTOXILYZER™ INSTRUMENT PRINTER CARD
© 1986 by CMI INC.

406
IS

SN 66-004798 E767.14 INVALID TEST	11/23/85 07:50
SAMPLE INTRODUCED AT IMPROPER TIME.	
<hr/> <hr/>	
SUBJECT'S NAME	
TIME FIRST OBSERVED	INSTRUMENT LOCATION
46331	
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	
<hr/> <hr/> <hr/>	
INTOXILYZER® INSTRUMENT PRINTER CARD	
© 1986 by CMI INC.	
	

586
IJ

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE

Agency Tucson Police Department Date 1/12/05 Time 0730
Intoxilyzer Serial # 66-004798 Location ODD
QA Specialist Irene Johnson 46331
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made.
Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓) 6. Push Start Test button or press enter on keyboard.
- (✓) 7. Insert card in response to display.
- (✓) 8. Air blank completed.
- (✓) 9. Calibration check completed.
Test results 0.078 & 0.079 AC
- (✓) 10. Air blank completed.
- (✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Irene D. Johnson 46331

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

TUCSON POLICE DEPARTMENT
INTOXICOLYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004798
11/23/05 07:52

TEST BRAC VALUE TIME
AIR BLANK .000 07:51
INVALID SAMPLE MOA .000 07:51
AIR BLANK .000 07:51

BEGIN NEW DEPRIVATION PERIOD

TUCSON POLICE DEPARTMENT
INTOXICOLYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004798
11/23/05

TEST BRAC VALUE TIME
AIR BLANK .000 07:51
AIR BLANK .000 07:51
*SUBJECT TEST :000 07:51
AIR BLANK .000 07:51

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

696
IJ

TEST BRAC VALUE TIME
AIR BLANK .000 07:51
AIR BLANK .000 07:51
*SUBJECT TEST :000 07:51
AIR BLANK .000 07:51

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME	ODD	INSTRUMENT LOCATION
TIME FIRST OBSERVED	46331	OPERATOR
ADDITIONAL INFORMATION AND/OR REMARKS		
Sweet breath		
INTOXICOLYZER® INSTRUMENT PRINTER CARD		
© 1986 by CMI INC.		

CMI Inc.