



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

je

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-004798 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 2.15.06 Cal 1 0.079 AC Cal 2 0.078 AC

Date: 3.16.06 Cal 1 0.078 AC Cal 2 0.079 AC

Date: _____ Cal 1 _____ AC Cal 2 _____ AC

Between the dates of 2.15.06 and 3.16.06
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.


Terry S. Gallegos

Terry S. Gallegos

3.23.2006

ADPS QAS Permit # QAS 976

SUBSCRIBED AND SWORN to before me this 23, day of
March, 2006.

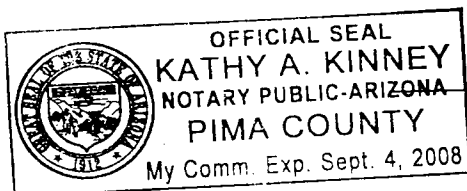

Kathy A. Kinney

Notary Public

*Includes documents specified in A.A.C. R9-14-404.A.6

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT P

106
754

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

gle

STANDARD CALIBRATION CHECK PROCEDURE

Agency Tucson Police Department Date 3.16.2006 Time 1127
Intoxilyzer Serial # 66-004798 Location ODD
QA Specialist Terry Gallegos 28847
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓) 6. Push Start Test button or press enter on keyboard.
- (✓) 7. Insert card in response to display.
- (✓) 8. Air blank completed.
- (✓) 9. Calibration check completed. Test results 0.078AC and 0.079 AC
- (✓) 10. Air blank completed.
- (✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Terry Gallegos

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

TEST	BRAC VALUE	TIME
AIR BLANK	.000	11:26
SUBJECT TEST	.000	11:27
AIR BLANK	.000	11:27

$$DVM = 319$$

SUBJECT'S NAME _____
TIME FIRST OBSERVED _____
INSTRUMENT LOCATION _____
OBD _____
Gallejos 28847 _____
APPLICATOR _____
LOCATION AND INFORMATION AND/OR REMARKS _____



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

[illegible]

SUBJECT'S NAME	INSTRUMENT LOCATION	TIME FIRST OBSERVED	OPERATOR	COORDINATE AND/OR REMARKS
	ODD		Gallegos	28847



© 1986 by CMI INC.

file 30f6 B4

Tucson Police Department Crime Laboratory

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
A2 MODEL 5000 SN 66-004798
03/16/06 11:31

TEST	BRAC VALUE	TIME
AIR BLANK	.000	11:30
CAL. CHECK	.079	11:30
AIR BLANK	.000	11:31

SUBJECT'S NAME
INSTRUMENT LOCATION
TIME FIRST OBSERVED

Gallegos 28847

OPERATOR
ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
A2 MODEL 5000 SN 66-004798
03/16/06 11:29

TEST	BRAC VALUE	TIME
AIR BLANK	.000	11:28
CAL. CHECK	.078	11:29
AIR BLANK	.000	11:29

SUBJECT'S NAME
INSTRUMENT LOCATION
TIME FIRST OBSERVED

Gallegos 28847

OPERATOR
ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC

4066
D4

SN 60-894798
E767.14
INVALID TEST

03/16/86
11:37

SAMPLE INTRODUCED
AT IMPROPER TIME.

SUBJECT'S NAME	
TIME FIRST OBSERVED	INSTRUMENT LOCATION
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	

INTOXILYZER* INSTRUMENT PRINTER CARD

© 1986 by CMI INC

CMI

EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

5066

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*AGENCY Tucson Police QA SPECIALISTTerry Gallegos 28847
(Print Name)

INTOXILYZER SERIAL #

66-004798

LOCATION

0DD

DATE

March 16182006

TIME

1127

☒ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

☒ DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .319.

☒ Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".

☒ Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".

☒ Clock time check. Mode selection switch S10 on or keyboard menu selection "E".

☒ Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

☒ Alcohol free subject Test result 0.000 AC

☒ Error recognition logic system functioning

Invalid test printed

☒ Proper sample recognition system

Invalid sample printed ☒

Deficient sample printed ☒

N/A Completeness of sample purge with collection tube**

☒ Calibration standard 0.080 AC

Results: 0.078 AC and 0.079 AC

Instrument operating properly and accurately. YES ☒ NO ☐

COMMENTS:

SIGNATURE

Terry Gallegos

0. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD
**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR
SAMPLE CAPTURE.

DHS/DSLS/Form C117 (Rev. 7-93)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

TUCSON POLICE DEPARTMENT
 INTOXILYZER - ALCOHOL ANALYZER
 A2 MODEL 5000 SN 66-004798
 03/16/06

TEST	BRAC VALUE	TIME
AIR BLANK	.000	11:37
#SUBJECT TEST	.225	11:40
AIR BLANK	.000	11:41

* DEFICIENT SAMPLE - VALUE
 PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME DDD

INSTRUMENT LOCATION

TIME FIRST OBSERVED Gallegos 28847

OPERATOR
 ADDITIONAL INFORMATION AND/OR REMARKS

Sweet breath



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
 INTOXILYZER - ALCOHOL ANALYZER
 A2 MODEL 5000 SN 66-004798
 03/16/06

TEST	BRAC VALUE	TIME
AIR BLANK	.000	11:35
INVALID SAMPLE	.000	11:36
AIR BLANK	.000	11:36

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME

DDD

INSTRUMENT LOCATION

TIME FIRST OBSERVED

Gallegos 28847

OPERATOR
 ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

60f6 TS4