



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

PS

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-004814 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 6/7/05 Cal 1 0.076 AC Cal 2 0.075 AC

Date: 7/7/05 Cal 1 0.077 AC Cal 2 0.077 AC

Date: — Cal 1 — AC Cal 2 — AC

Between the dates of 6/7/05 and 7/7/05 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on NA.
The instrument was returned to service on NA.

Irene Doran

Irene Doran

7/19/05

ADPS QAS Permit # QAS 29586
SUBSCRIBED AND SWORN to before me this 19, day of
July, 2005.

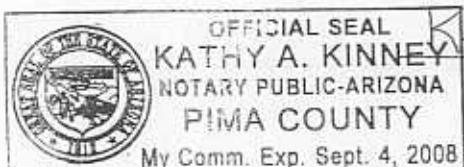

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R9-14-404.A.6

1076
P4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000* P4

AGENCY Tucson Police QA SPECIALIST Trene Doran 46331
 (Print Name)
 INTOXILYZER SERIAL # 66-004814 LOCATION ODM
 DATE July 7 2005 TIME 1231

Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .271.
- Display Test check. Mode selection switch S1 on and S2 and S3 off or Keyboard menu selection "V".
- Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- Alcohol free subject Test result 0.000 AC
 - Error recognition logic system functioning
 - Invalid test printed
 - Proper sample recognition system
 - Invalid sample printed
 - Deficient sample printed
 - N/A Completeness of sample purge with collection tube**
 - Calibration standard 0.080 AC
 - Results: 0.077AC + 0.077AC
- Instrument operating properly and accurately. YES NO

COMMENTS: * Moved clock ahead 2 min.

SIGNATURE: Trene Doran

0. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR SAMPLE CAPTURE.

SN 66-004814	07/07/05
E767.14	07/12:34
EEEEE EEEEE EEEEE EEEEE	
HHHHH HHHHH HHHHH HHHHH	
0307	-----
0306	-----
0305	-----
0302	-----
0303	-----
0304	-----
0301	-----

()

TUCSON POLICE DEPARTMENT	INTOXILYZER -	ALCOHOL ANALYZER
INTOXILYZER -	SH 66-004814	SH 66-004814
AZ MODEL 5000	12:33	
07/07/05		
TEST	BRAC VALUE	TIME
AIR BLANK	.000	12:32
SUBJECT TEST	.000	12:32
AIR BLANK	.000	12:33

DVM = .271

SUBJECT'S NAME	ODM
INSTRUMENT LOCATION	
TIME FIRST OBSERVED	
OPERATOR	46331
ADDITIONAL INFORMATION AND/OR REMARKS	

286
6/1

INSTRUMENT PRINTER CARD	INTOXILYZER® INSTRUMENT PRINTER CARD
© 1986 by CMI INC.	© 1986 by CMI INC.

CMI

386
P4

Tucson Police Department Crime Laboratory

TUCSON POLICE DEPARTMENT		DEPARTMENT	
INTOXILYZER - ALCOHOL ANALYZER		INTOXILYZER - ALCOHOL ANALYZER	
AZ MODEL 5000	SN 66-004814	AZ MODEL 5000	SN 66-004814
07/07/05	12:46	07/07/05	12:47
TEST	BRAC VALUE	TEST	BRAC VALUE
AIR BLANK	.000	AIR BLANK	.000
CAL. CHECK	.077	CAL. CHECK	.077
AIR BLANK	.000	AIR BLANK	.000
	12:45		12:46
	12:45		12:47
	12:45		12:47

SUBJECT'S NAME	ODM
INSTRUMENT LOCATION	OFFICE
TIME FIRST OBSERVED	46331
OPERATOR	ODM
ADDITIONAL INFORMATION AND/OR REMARKS	

INTOXILYZER® INSTRUMENT PRINTER CARD	© 1986 by CMI INC.
CMI	

4076
PA

SN 66-004814 E767.14 INVALID TEST	07/07/05 12:44
SAMPLE INTRODUCED AT IMPROPER TIME.	
<hr/> <hr/>	
SUBJECT'S NAME	
TIME FIRST OBSERVED	
INSTRUMENT LOCATION	
40331	
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	
<hr/> <hr/> <hr/>	
INTOXILYZER® INSTRUMENT PRINTER CARD	
© 1986 by CMI INC.	
CMI	

596
45

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*STANDARD CALIBRATION CHECK PROCEDURE *BH*

Agency Tucson Police Department Date 7/7/05 Time 1231
Intoxilyzer Serial # 66-004814 Location ODM
QA Specialist Irene Doran 40331
(Print Name)

- (1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made.
Standard value: 0.080 AC
- (2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$.
- (3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (6. Push Start Test button or press enter on keyboard.
- (7. Insert card in response to display.
- (8. Air blank completed.
- (9. Calibration check completed.
Test results 0.077 + 0.077 AC
- (10. Air blank completed.
- (11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Irene Doran

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004814
07/07/95 12:34

TEST	BRAC VALUE	TIME
AIR BLANK	.000	12:33
INVALID SAMPLE	:NOA	12:34
AIR BLANK	.000	12:34

BEGIN NEW DEPRIVATION PERIOD

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004814
07/07/95 12:34

TEST	BRAC VALUE	TIME
AIR BLANK	.000	12:33
*SUBJECT TEST		
AIR BLANK	.000	12:43

TEST
AIR BLANK : .000 12:33
*SUBJECT TEST : .000 12:43
AIR BLANK : .000 12:43

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

696
25

SUBJECT'S NAME	S. B. ECT'S NAME
TIME FIRST OBSERVED	TIME FIRST OBSERVED
INSTRUMENT LOCATION	INSTRUMENT LOCATION
OPERATOR	OPERATOR
ADDITIONAL INFORMATION AND/OR REMARKS	
INTOXILYZER® INSTRUMENT PRINTER CARD	
© 1986 by CMI INC.	

46331
Sweet breath

CMI