



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

Handwritten signature

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-004814 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 5.11.06 Cal 1 0.077 AC Cal 2 0.077 AC

Date: 6.8.06 Cal 1 0.081 AC Cal 2 0.080 AC

Date: — Cal 1 — AC Cal 2 — AC

Between the dates of 5.11.06 and 6.8.06 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Handwritten signature of Terry S. Gallegos
(Terry S. Gallegos)

6.13.2006

ADPS QAS Permit # QAS 976

SUBSCRIBED AND SWORN to before me this 13, day of June, 2006.

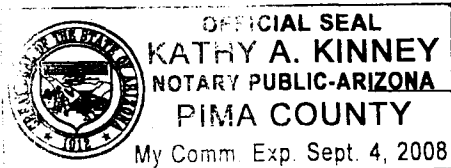

Handwritten signature of Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R9-14-404.A.6

Note: The applicable A.A.C. document sections were revised effective May 18, 2006 at 0900 hours. The new sections are A.A.C. R13-10-104.A.5 and A.A.C. R13-10-104.A.2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1056
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry S. Gallegos 28847 AGENCY Tucson Police Department
DATE 06-08-2006 TIME 1402
INTOXILYZER SERIAL # 66-004814 LOCATION ODM

- (✓) 1. Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- (✓) 1. DVM Test check. Setting should be between .010 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is 0.303.
- (✓) 2. Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- (✓) 3. Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- * (✓) 4. Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- (✓) 5. Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject Test result 0. 000 AC
- (✓) 2. Error recognition logic system functioning
Invalid test printed
- (✓) 3. Proper sample recognition system
Invalid sample printed ✓
Deficient sample printed ✓
- (✓) 4. Standard Calibration Check standard 0. 080 AC
Results: 0. 081AC and 0.080 AC

Instrument operating accurately and properly. YES ✓ NO

COMMENTS: * 1354 -> 1401

SIGNATURE Terry S. Gallegos

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DPS Form Exh E-4 (Rev. 05-1)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

$$DVM = 0.303$$

SUBJECT'S NAME _____
TIME FIRST OBSERVED _____
INSTRUMENT LOCATION _____
ODM _____
OPERATOR _____
ADDITIONAL INFORMATION AND/OR REMARKS _____

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI INC.

2086 TBH Jle

SUBJECT'S NAME	TIME FIRST OBSERVED	INSTRUMENT LOCATION	OPERATOR	ADDITIONAL INFORMATION AND/OR REMARKS
		DDM	Gallegos	28847

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI INC.

Tucson Police Department Crime Laboratory

3066
TS4 file

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004814
06/08/06 14:06

TEST BRAC VALUE TIME
AIR BLANK .000 14:05
CAL. CHECK .000 14:06
AIR BLANK .000 14:06

SUBJECT'S NAME
TIME FIRST OBSERVED
INSTRUMENT LOCATION

Gallegos 28847

OPERATOR
ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004814
06/08/06 14:04

TEST BRAC VALUE TIME
AIR BLANK .000 14:04
CAL. CHECK .001 14:04
AIR BLANK .000 14:04

SUBJECT'S NAME
TIME FIRST OBSERVED
INSTRUMENT LOCATION

Gallegos 28847

OPERATOR
ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

4 of 6
TS4

file

SN 66-004814
E767.14
INVALID TEST

06/08/06
14:17

SAMPLE INTRODUCED
AT IMPROPER TIME.

SUBJECT'S NAME	
ODM	
TIME FIRST OBSERVED	INSTRUMENT LOCATION
Gallegos	28847
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	

INTOXILYZER[®] INSTRUMENT PRINTER CARD

© 1986 by CMI INC.



THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

506
Dh

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST Terry S. Gallegos 28847 AGENCY Tucson Police Department
DATE 06-08-2006 TIME 1402
INTOXILYZER SERIAL # 66-004814 LOCATION DDM

- (✓) 1. Pour a standard alcohol concentration solution into a clean dry simulator jar and assemble the simulator. Ensure that a tight seal has been made. Standard value: 0. DBO AC.
- (✓) 2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓✓) 6. Push Start Test button or press enter on keyboard.
- (✓✓) 7. Insert card in response to display.
- (✓✓) 8. Air blank completed.
- (✓✓) 9. Standard Calibration Check completed. Test results 0. DB1 AC and 0. DBO AC.
- (✓✓) 10. Air blank completed.
- (✓✓) 11. When test complete, remove printed record. Attach the record to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

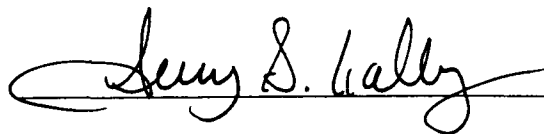

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DPS Form Exh E-2 (Rev. 05-1)

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004814
06/08/06 14:16

TEST	BRAC VALUE	TIME
AIR BLANK	.000	14:13
*SUBJECT TEST	.049	14:16
AIR BLANK	.000	14:16

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME

ODM

INSTRUMENT LOCATION

TIME FIRST OBSERVED

Gallegos 20847

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

CMI

INTOXILYZER* INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004814
06/08/06 14:12

TEST	BRAC VALUE	TIME
AIR BLANK	.000	14:12
INVALID SAMPLE	.MOA	14:12
AIR BLANK	.000	14:12

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME

ODM

INSTRUMENT LOCATION

TIME FIRST OBSERVED

Gallegos 20847

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet breath

CMI

INTOXILYZER* INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

6066
TS4
je