



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

fm

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-064922 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 11/23/05 Cal 1 0.076 AC Cal 2 0.077 AC

Date: 12/20/05 Cal 1 0.077 AC Cal 2 0.077 AC

Date: - Cal 1 - AC Cal 2 - AC

Between the dates of 11/23/05 and 12/20/05 the instrument was working properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Irene Doran Johnson

Irene Doran Johnson

12/22/05

ADPS QAS Permit # QAS 29586

SUBSCRIBED AND SWORN to before me this 22, day of December, 2005.

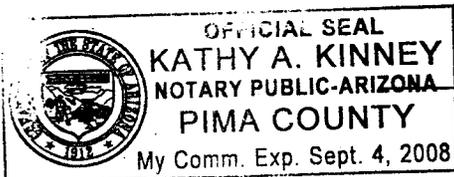

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R9-14-404.A.6

EXHIBIT Q

196
EJ

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES *BN*
INTOXILYZER MODEL 5000*

AGENCY Tucson Police QA SPECIALIST Irene Johnson 40331
(Print Name)

INTOXILYZER SERIAL # 60-004922 LOCATION ODE
DATE December 20 18 2005 TIME 1155

Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .347.
- Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- Alcohol free subject Test result 0.000 AC
 - Error recognition logic system functioning
Invalid test printed
 - Proper sample recognition system
Invalid sample printed
Deficient sample printed
 - N/A Completeness of sample purge with collection tube**
 - Calibration standard 0.080 AC
Results: 0.077AC + 0.077AC
- Instrument operating properly and accurately. YES NO

COMMENTS: _____

SIGNATURE Irene Johnson 40331

0. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD
**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR SAMPLE CAPTURE.

306
ID

Tucson Police Department Crime Laboratory

TUCSON POLICE DEPARTMENT
 INTOXILYZER - ALCOHOL ANALYZER
 AZ MODEL 5000 SN 66-004922
 12/20/05 12:08

TEST BRAC VALUE TIME
 AIR BLANK .000 12:07
 CAL. CHECK .077 12:07
 AIR BLANK .000 12:08

SUBJECT'S NAME

ODE

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46331

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
 INTOXILYZER - ALCOHOL ANALYZER
 AZ MODEL 5000 SN 66-004922
 12/20/05 12:06

TEST BRAC VALUE TIME
 AIR BLANK .000 12:05
 CAL. CHECK .077 12:06
 AIR BLANK .000 12:06

SUBJECT'S NAME

ODE

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46331

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

406
EJ

SN 66-004922
E767.14
INVALID TEST

12/20/05
12:01

SAMPLE INTRODUCED
AT IMPROPER TIME.



SUBJECT'S NAME

ODE

TIME FIRST OBSERVED

INSTRUMENT LOCATION

40331

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER® INSTRUMENT PRINTER CARD



© 1986 by CMI INC.

EXHIBIT P

566
EJ

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE *AM*

Agency Tucson Police Department Date 12/20/05 Time 1155
Intoxilyzer Serial # 66-004922 Location ODE
QA Specialist Irene Johnson 46331
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach 34°C+.2°C.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓) 6. Push Start Test button or press enter on keyboard.
- (✓) 7. Insert card in response to display.
- (✓) 8. Air blank completed.
- (✓) 9. Calibration check completed. Test results 0.077+0.077 AC
- (✓) 10. Air blank completed.
- (✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE Irene Johnson 46331

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004922
12/20/05

TEST BRAC VALUE TIME
AIR BLANK .000 11:59
INVALID SAMPLE .MOA 12:00
AIR BLANK .000 12:00

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME

ODE

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46331

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet breath



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004922
12/20/05

TEST BRAC VALUE TIME
AIR BLANK .000 12:01
*SUBJECT TEST .200 12:04
AIR BLANK .000 12:04

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME

ODE

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46331

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

986
FI