



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

file

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-004922 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 2·15·06 Cal 1 0.081 AC Cal 2 0.080 AC

Date: 3·16·06 Cal 1 0.075 AC Cal 2 0.077 AC

Date: _____ Cal 1 _____ AC Cal 2 _____ AC

Between the dates of 2·15·06 and 3·16·06, the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

Terry S. Gallegos

Terry S. Gallegos


3·23·2006

ADPS QAS Permit # QAS 976

SUBSCRIBED AND SWORN to before me this 23, day of
March, 2006.

Kathy A Kinney

Notary Public

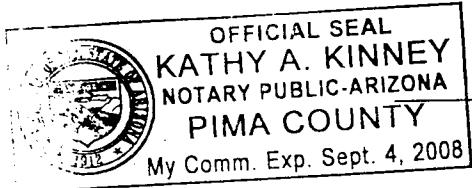

*Includes documents specified in A.A.C. R9-14-404.A.6

EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404 (A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*10ft
T84

AGENCY Tucson Police QA SPECIALIST Terry Gallegos 28947
(Print Name)

INTOXILYZER SERIAL # 66-004922 LOCATION ODE
DATE March 16 TIME 1829

Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .368.
- Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- Alcohol free subject Test result 0.000 AC
 - Error recognition logic system functioning
 - Invalid test printed
 - Proper sample recognition system
 - Invalid sample printed
 - Deficient sample printed
 - N/A Completeness of sample purge with collection tube**
 - Calibration standard 0.080 AC
 - Results: 0.075AC and 0.077AC
- Instrument operating properly and accurately. YES NO

COMMENTS:

SIGNATURE Terry Gallegos

0. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR
SAMPLE CAPTURE.

SN 66-004922 DEPARTMENT
 INTOXILYZER - ALCOHOL ANALYZER
 AZ MODEL 5000 SN 66-004922
 03/16/86 18:30

	TEST	BRAC VALUE	TIME
E767.14	AIR BLANK	.000	18:29
EEEEEEEEE	SUBJECT TEST	.000	18:29
EEEEEEEEE	AIR BLANK	.000	18:30
Q386	-	-	-
Q385	-	-	-
Q382	-	-	-
Q383	-	-	-
Q384	-	-	-
Q381	-	-	-

DIVM = .368

SUBJECT'S NAME: ODE
 INSTRUMENT LOCATION: OPERATOR
 TIME FIRST OBSERVED: ODE
 ADDITIONAL INFORMATION AND/OR REMARKS: Gall ego 5 28847

SUBJECT'S NAME: ODE
 INSTRUMENT LOCATION: OPERATOR
 TIME FIRST OBSERVED: ODE
 ADDITIONAL INFORMATION AND/OR REMARKS: Gall ego 5 28847

INTOXILYZER INSTRUMENT PRINTER CARD
 © 1986 by CMI INC.

CMI Inc.

INTOXILYZER INSTRUMENT PRINTER CARD
 © 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004922
03/16/86 18:41

TEST	BRAC VALUE	TIME
AIR BLANK	:000	18:40
CAL. CHECK	:075	18:41
AIR BLANK	:000	18:41

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004922
03/16/86 18:43

TEST	BRAC VALUE	TIME
AIR BLANK	:000	18:40
CAL. CHECK	:075	18:41
AIR BLANK	:000	18:43

Tucson Police Department Crime Laboratory

30f6
PS4
JL

SUBJECT'S NAME	ODE
INSTRUMENT LOCATION	DRATON
TIME FIRST OBSERVED	
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	
Gallegos 20847	

SUBJECT'S NAME	ODE
INSTRUMENT LOCATION	DRATON
TIME FIRST OBSERVED	
OPERATOR	
ADDITIONAL INFORMATION AND/OR REMARKS	
Gallegos 20847	

INTOXILYZER INSTRUMENT PRINTER CARD
© 1986 by CMI INC.

CMI INC

INTOXILYZER INSTRUMENT PRINTER CARD
© 1986 by CMI INC.

CMI INC

Hofg
TSG

SN 66-004922 03/16/86
E767.14 18:45
INVALID TEST

TEST	BAC VALUE	TIME
AIR BLANK	.000	18:44
SUBJECT TEST	.000	18:45

SAMPLE INTRODUCED
AT IMPROPER TIME.

SUBJECT'S NAME

ODE

TIME FIRST OBSERVED

INSTRUMENT LOCATION

Gullegos 28847

OPERATOR

ADDITIONAL INFORMATION AND / OR REMARKS

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

5056
TS4

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE

Agency Tucson Police Department Date 3/16/2006 Time 1829
 Intoxilyzer Serial # 66-004922 Location ODE
 QA Specialist Terry Gallegos 28847
 (Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator.
Insure that a tight seal has been made.
Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓✓) 6. Push Start Test button or press enter on keyboard.
- (✓✓) 7. Insert card in response to display.
- (✓✓) 8. Air blank completed.
- (✓✓) 9. Calibration check completed.
Test results 0.075AC and 0.077 AC
- (✓✓) 10. Air blank completed.
- (✓✓) 11. When display reads Test Complete remove evidence card.
Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Jerry Gallegos

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004922
03/16/06

TEST BRAC VALUE TIME
AIR BLANK .000 18:34
INVALID SAMPLE : NOA 18:34
AIR BLANK .000 18:34

BEGIN NEW DEPRIVATION PERIOD

TUCSON POLICE DEPARTMENT DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-004922
03/16/06

TEST BRAC VALUE TIME
AIR BLANK .000 18:34
*SUBJECT TEST : .005 18:38
AIR BLANK .000 18:38

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME	ODE
INSTRUMENT LOCATION	TIME FIRST OBSERVED
OPERATOR	ADDITIONAL INFORMATION AND/OR REMARKS

ODE

TIME FIRST OBSERVED

ADDITIONAL INFORMATION AND/OR REMARKS

SUBJECT'S NAME	ODE
INSTRUMENT LOCATION	TIME FIRST OBSERVED
OPERATOR	ADDITIONAL INFORMATION AND/OR REMARKS

ODE

TIME FIRST OBSERVED

ADDITIONAL INFORMATION AND/OR REMARKS

SUBJECT'S NAME	ODE
INSTRUMENT LOCATION	TIME FIRST OBSERVED
OPERATOR	ADDITIONAL INFORMATION AND/OR REMARKS

ODE

TIME FIRST OBSERVED

ADDITIONAL INFORMATION AND/OR REMARKS

SUBJECT'S NAME	ODE
INSTRUMENT LOCATION	TIME FIRST OBSERVED
OPERATOR	ADDITIONAL INFORMATION AND/OR REMARKS

ODE

TIME FIRST OBSERVED

ADDITIONAL INFORMATION AND/OR REMARKS

b6f6
TSY

© 1986 by CMI INC.



INTOXILYZER® INSTRUMENT PRINTER CARD



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.