



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

154

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 5000 Serial Number 66-005003 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 7/13/05 Cal 1 0.079 AC Cal 2 0.078 AC

Date: 8/10/05 Cal 1 0.079 AC Cal 2 0.079 AC

Date: — Cal 1 — AC Cal 2 — AC

Between the dates of 7/13/05 and 8/10/05 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


Irene Doran

Irene Doran

8/15/05

ADPS QAS Permit # QAS 29586

SUBSCRIBED AND SWORN to before me this 15, day of August, 2005.

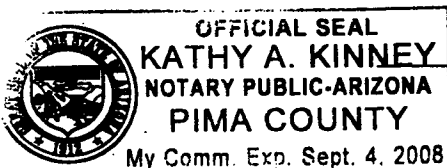

Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R9-14-404.A.6

EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 5000\*

AGENCY Tucson Police QA SPECIALIST Irene Doran 46331  
(Print Name)

INTOXILYZER SERIAL # 66-005003 LOCATION DUI HQ  
DATE August 10 19 2005 TIME 1036

☒ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

☒ DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is 0.351.

☒ Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".

☒ Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".

☒ Clock time check. Mode selection switch S10 on or keyboard menu selection "E".

☒ Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

☒ Alcohol free subject Test result 0.000 AC

☒ Error recognition logic system functioning  
Invalid test printed

☒ Proper sample recognition system  
Invalid sample printed ☒  
Deficient sample printed ☒

N/A Completeness of sample purge with collection tube\*\*

☒ Calibration standard 0.080 AC

Results: 0.079AC + 0.079AC

Instrument operating properly and accurately. YES ☒ NO ☐

COMMENTS:

SIGNATURE

Irene Doran  
0. N/A alcohol concentration standard collected for subsequent analysis.\*\*

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD  
\*\*THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR  
SAMPLE CAPTURE.

DHS/DSLS/Form C117 (Rev. 7-93)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



306  
ps

# Tucson Police Department Crime Laboratory

CMI INC  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-005003  
08/10/05 10:55

TEST	BRAC VALUE	TIME
AIR BLANK	.000	10:54
CAL. CHECK	.079	10:54
AIR BLANK	.000	10:55

SUBJECT'S NAME

DUI #0  
INSTRUMENT LOCATION

TIME FIRST OBSERVED

40331

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI INC  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-005003  
08/10/05 10:53

TEST	BRAC VALUE	TIME
AIR BLANK	.000	10:52
CAL. CHECK	.079	10:53
AIR BLANK	.000	10:53

SUBJECT'S NAME

DUI #0  
INSTRUMENT LOCATION

TIME FIRST OBSERVED

40331

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

4086  
PE

SN 66-005303  
E767.14  
INVALID TEST

08/10/05  
10:47

SAMPLE INTRODUCED  
AT IMPROPER TIME.

SUBJECT'S NAME

Dut HQ

TIME FIRST OBSERVED

INSTRUMENT LOCATION

46331

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER<sup>®</sup> INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

**CMI**  
INC.

EXHIBIT P

506  
PS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 5000\*

STANDARD CALIBRATION CHECK PROCEDURE

AM

Agency Tucson Police Department Date 8/10/05 Time 1036  
Intoxilyzer Serial # 66-005003 Location DUI HQ  
QA Specialist Irene Doran 46331  
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach  $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$ .
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓) 6. Push Start Test button or press enter on keyboard.
- (✓) 7. Insert card in response to display.
- (✓) 8. Air blank completed.
- (✓) 9. Calibration check completed. Test results 0.079 AC + 0.079 AC
- (✓) 10. Air blank completed.
- (✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Irene Doran

\*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

CMI INC  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-005003  
08/10/05 10:42

TEST BRAC VALUE TIME  
AIR BLANK .000 10:41  
INVALID SAMPLE .000 10:42  
AIR BLANK .000 10:42

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME

DUI HQ  
INSTRUMENT LOCATION

TIME FIRST OBSERVED

40331  
OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet breath



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI INC  
INTOXILYZER - ALCOHOL ANALYZER  
AZ MODEL 5000 SN 66-005003  
08/10/05 10:46

TEST BRAC VALUE TIME  
AIR BLANK .000 10:43  
\*SUBJECT TEST .123 10:46  
AIR BLANK .000 10:46

\* DEFICIENT SAMPLE - VALUE  
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME

DUI HQ  
INSTRUMENT LOCATION

TIME FIRST OBSERVED

40331  
OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

696  
PK