



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

BM

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-005004 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R9-14-404.A.3, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 9/8/05 Cal 1 0.081 AC Cal 2 0.080 AC

Date: 10/6/05 Cal 1 0.079 AC Cal 2 0.080 AC

Date: N/A Cal 1 AC Cal 2 AC

Between the dates of 9/8/05 and 10/6/05 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.


The instrument was returned to service on N/A.

Jennifer L. Coombs
Jennifer L. Coombs

10/14/05

ADPS QAS Permit # QAS 29886

SUBSCRIBED AND SWORN to before me this 14, day of October, 20 05.

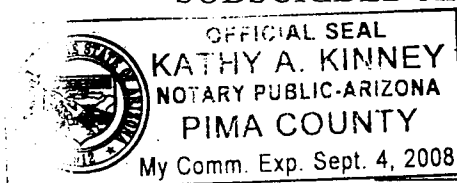

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R9-14-404.A.6

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT Q

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

AGENCY Tucson Police QA SPECIALIST Jennifer Combs 46160
(Print Name)

INTOXILYZER SERIAL # 66-005204 LOCATION COOMBS' HOME
DATE October 6 18 2005 TIME 1604

☒ Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

☒ DVM Test check. Setting should be between .200 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .386.

☒ Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".

☒ Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".

☒ Clock time check. Mode selection switch S10 on or keyboard menu selection "E".

☒ Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

☒ Alcohol free subject Test result 0.000 AC

☒ Error recognition logic system functioning
Invalid test printed

☒ Proper sample recognition system
Invalid sample printed ✓
Deficient sample printed ✓

N/A Completeness of sample purge with collection tube**

☒ Calibration standard 0.080 AC
Results: 0.079 AC & 0.080 AC

Instrument operating properly and accurately. YES ☒ NO ☐

COMMENTS: *The instrument was dropped off at DUI HQ as I was leaving. I brought the instrument home and completed the SDAP there so I could return it to the DUI HQ the next morning.

SIGNATURE Jennifer L. Combs

0. N/A alcohol concentration standard collected for subsequent analysis.**

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD
**THIS STEP IS ONLY REQUIRED IF THIS DEVICE IS BEING USED FOR
SAMPLE CAPTURE.

DHS/DSLS/Form C117 (Rev. 7-93)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

1 of 6
fle

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:11
SUBJECT TEST	.000	16:12
AIR BLANK	.000	16:13

SUBJECT'S NAME _____

INSTRUMENT LOCATION

46160

ADDITIONAL INFORMATION AND/OR REMARKS



INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

2 of 6
file

SUBJECT'S NAME _____

INSTRUMENT LOCATION

4/6/60

[illegible]

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

Tucson Police Department Crime Laboratory *py*

3 of 6
fle

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-005004
10/06/05 16:07

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:09
CAL. CHECK	.000	16:10
AIR BLANK	.000	16:10

SUBJECT'S NAME

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160
OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER™ INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-005004
10/06/05 16:07

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:06
CAL. CHECK	.079	16:07
AIR BLANK	.000	16:07

SUBJECT'S NAME

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160
OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER™ INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

4 of 6
file

SN 66-035334
E767.14
INVALID TEST

10/26/85
16:13

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:13

SAMPLE INTRODUCED
AT IMPROPER TIME.

SUBJECT'S NAME

TIME FIRST OBSERVED

INSTRUMENT LOCATION

46160
OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER[®] INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI^{INC.}

EXHIBIT P

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R9-14-404(A)

ARIZONA DEPARTMENT OF HEALTH SERVICES
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE *BH*

Agency Tucson Police Department Date 10/6/05 Time 1604
Intoxilyzer Serial # 66-005004 Location COOMBS' HOME
QA Specialist Jennifer Coombs 46160
(Print Name)

- (✓) 1. Pour a standard alcohol solution of known value into a clean dry simulator jar and assemble the simulator. Insure that a tight seal has been made. Standard value: 0.080 AC
- (✓) 2. Turn on the simulator and allow the temperature to reach $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- (✓) 6. Push Start Test button or press enter on keyboard.
- (✓) 7. Insert card in response to display.
- (✓) 8. Air blank completed.
- (✓) 9. Calibration check completed. Test results 0.079AC # 0.080 AC
- (✓) 10. Air blank completed.
- (✓) 11. When display reads Test Complete remove evidence card. Attach the card to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE

Jennifer L. Coombs

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DHS/DSLS/Form C116 (Rev. 7-93)

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-005084
10/06/05 16:14

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:14
INVALID SAMPLE	.000	16:14
AIR BLANK	.000	16:14

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet Breath

INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.



TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-005031
10/06/05 16:19

TEST	BRAC VALUE	TIME
AIR BLANK	.000	16:19
*SUBJECT TEST	.119	16:19
AIR BLANK	.000	16:19

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER INSTRUMENT PRINTER CARD

© 1986 by CMI INC.



6 of 6
file