



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

184

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 166-005004 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 7/13/06 Cal 1 0.079 AC Cal 2 0.081 AC

Date: 8/12/06 Cal 1 0.079 AC Cal 2 0.080 AC

Date: N/A Cal 1 _____ AC Cal 2 _____ AC


Between the dates of 7/13/06 and 8/12/06 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
 The instrument was returned to service on N/A.

Jennifer L. Coombs
 Jennifer L. Coombs

8/22/06

ADPS QAS Permit # QAS 29886
 SUBSCRIBED AND SWORN to before me this 22, day of August, 2006.

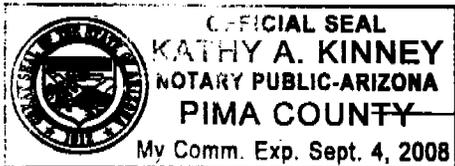

Kathy A. Kinney
 Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD QUALITY ASSURANCE PROCEDURE

1 off file
BH

QA SPECIALIST Jennifer Coombs 46160 AGENCY Tucson Police Department
DATE August 12, 2006 TIME 1140
INTOXILYZER SERIAL # 106-005004 LOCATION DUI HQ

- () 1. Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- () 1. DVM Test check. Setting should be between .010 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .376.
- () 2. Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- () 3. Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- () 4. Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- () 5. Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- () 1. Alcohol-free subject Test result 0. 000 AC
- () 2. Error recognition logic system functioning
Invalid test printed
- () 3. Proper sample recognition system
Invalid sample printed ✓
Deficient sample printed ✓
- () 4. Standard Calibration Check standard 0. 080 AC
Results: 0. 079AC & 0.080 AC

Instrument operating accurately and properly. YES NO

COMMENTS:

SIGNATURE Jennifer L. Coombs

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

Tucson Police Department Crime Laboratory

BM
3 of 6
file

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-0055004
08/12/86 11:51

TEST BRAC VALUE TIME
AIR BLANK .000 11:50
CAL. CHECK .000 11:50
AIR BLANK .000 11:51

SUBJECT'S NAME
INSTRUMENT LOCATION
DUI HQ

TIME FIRST OBSERVED

OPERATOR
46160

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-0055004
08/12/86 11:49

TEST BRAC VALUE TIME
AIR BLANK .000 11:48
CAL. CHECK .079 11:49
AIR BLANK .000 11:49

SUBJECT'S NAME
INSTRUMENT LOCATION
DUI HQ

TIME FIRST OBSERVED

OPERATOR
46160

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

4 of 6
file

SN 66-005004
E767.14
INVALID TEST

08/13/01
11:45

PRESSED START TEST AT
IMPROPER TIME

SUBJECT'S NAME

DUI HQ

TIME FIRST OBSERVED

INSTRUMENT LOCATION

OPERATOR

46160

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER® INSTRUMENT PRINTER CARD



© 1986 by CMI INC.

EXHIBIT E-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

5 of 6
file

STANDARD CALIBRATION CHECK PROCEDURE

QA SPECIALIST Jennifer Coombs 4611aD AGENCY Tucson Police Department
DATE August 12 2006 TIME 1140
INTOXILYZER SERIAL # 66-005004 LOCATION DUI HQ

- (✓) 1. Pour a standard alcohol concentration solution into a clean dry simulator jar and assemble the simulator. Ensure that a tight seal has been made. Standard value: 0. 080 AC.
- (✓) 2. Turn on the simulator and allow the temperature to reach 34°C ± 0.2°C.
- (✓) 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- (✓) 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- ✓ (✓) 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- ✓ (✓) 6. Push Start Test button or press enter on keyboard.
- ✓ (✓) 7. Insert card in response to display.
- ✓ (✓) 8. Air blank completed.
- ✓ (✓) 9. Standard Calibration Check completed. Test results 0. 079 AC & 0. 080 AC.
- ✓ (✓) 10. Air blank completed.
- ✓ (✓) 11. When test complete, remove printed record. Attach the record to the completed checklist.
- (✓) 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE Jennifer L Coombs

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DPS Form Exh E-2 (Rev. 05-1)

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-005504
08/12/86 11:46

TEST BRAC VALUE TIME
AIR BLANK .000 11:45
INVALID SAMPLE .000 11:45
AIR BLANK .000 11:46

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME

DUI HQ

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet Breath

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI INC.

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-005504
08/12/86 11:55

TEST BRAC VALUE TIME
AIR BLANK .000 11:52
*SUBJECT TEST .147 11:54
AIR BLANK .000 11:55

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME

DUI HQ

INSTRUMENT LOCATION

TIME FIRST OBSERVED

46160

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

CMI INC.

73 606 fle