



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 5000 Serial Number 66-005004 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.080AC solution are 0.070AC to 0.090AC.

Date: 8.12.06 Cal 1 0.079 AC Cal 2 0.080 AC

Date: 9.8.06 Cal 1 0.079 AC Cal 2 0.079 AC

Date: N/A Cal 1 _____ AC Cal 2 _____ AC

Between the dates of 8.12.06 and 9.8.06 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.


Terry S. Gallegos

Terry S. Gallegos

9.12.2006

ADPS QAS Permit # QAS 976

SUBSCRIBED AND SWORN to before me this 12, day of September, 2006.

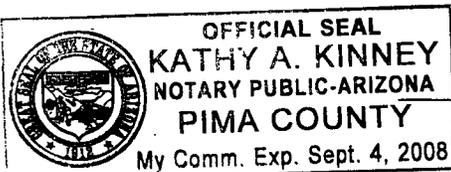

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1056
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 9-8-2006 TIME 1635
INTOXILYZER SERIAL # 66-005004 LOCATION DUI HQ

() 1. Display reads "READY TO START" or "PUSH BUTTON TO START TEST."

DIAGNOSTIC TESTS

- () 1. DVM Test check. Setting should be between .010 and .600. Mode selection switch S2 on, S1 and S3 off or keyboard menu selection "H". Reading is .365.
- () 2. Display Test check. Mode selection switch S1 on and S2 and S3 off or keyboard menu selection "V".
- () 3. Printer Test check. Mode selection switch S1, S2, S3 off or keyboard menu selection "P".
- () 4. Clock time check. Mode selection switch S10 on or keyboard menu selection "E".
- () 5. Date check. Mode selection switch S11 on or keyboard menu selection "E".

OPERATIONAL TESTS

- () 1. Alcohol-free subject Test result 0.000 AC
- () 2. Error recognition logic system functioning
Invalid test printed
- () 3. Proper sample recognition system
Invalid sample printed ✓
Deficient sample printed ✓
- () 4. Standard Calibration Check standard 0.080 AC
Results: 0.079AC and 0.079 AC

Instrument operating accurately and properly. YES NO

COMMENTS:

SIGNATURE

Terry Gallegos

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DPS Form Exh E-4 (Rev. 05-1)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

Tucson Police Department Crime Laboratory

3086
Tsh
JL

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-085834
09/08/06 16:36

TEST BRAC VALUE TIME
AIR BLANK .000 16:35
CAL. CHECK .079 16:35
AIR BLANK .000 16:36

SUBJECT'S NAME

DUI HQ

INSTRUMENT LOCATION

TIME FIRST OBSERVED

Gallegos

OPERATOR

20847

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.



TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-085834
09/08/06 16:37

TEST BRAC VALUE TIME
AIR BLANK .000 16:37
CAL. CHECK .079 16:37
AIR BLANK .000 16:37

SUBJECT'S NAME

DUI HQ

INSTRUMENT LOCATION

TIME FIRST OBSERVED

Gallegos

OPERATOR

20847

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.



4076 B4

je

SN 66-305884
E797.14
INVALID TEST

03/08/88
16:47

SAMPLE INTRODUCED
AT IMPROPER TIME.

SUBJECT'S NAME

DUI HQ

INSTRUMENT LOCATION

TIME FIRST OBSERVED

Gallegos 28847

OPERATOR

ADDITIONAL INFORMATION AND/OR REMARKS

INTOXILYZER[®] INSTRUMENT PRINTER CARD

CMI INC.

© 1986 by CMI INC.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

*Soft
TS4*

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 5000*

STANDARD CALIBRATION CHECK PROCEDURE

file

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 9.8.2006 TIME 1630
INTOXILYZER SERIAL # 66-005004 LOCATION DUI HQ

- () 1. Pour a standard alcohol concentration solution into a clean dry simulator jar and assemble the simulator. Ensure that a tight seal has been made. Standard value: 0.080 AC.
- () 2. Turn on the simulator and allow the temperature to reach 34°C ± 0.2°C.
- () 3. Set Intoxilyzer mode selection in the ACA mode by switching mode selection switch #9 on or selecting "C" on keyboard menu.
- () 4. Attach simulator to the simulator entrance port on the Intoxilyzer.
- () 5. Intoxilyzer 5000 display reads "READY TO START" or "PUSH BUTTON".
- () 6. Push Start Test button or press enter on keyboard.
- () 7. Insert card in response to display.
- () 8. Air blank completed.
- () 9. Standard Calibration Check completed. Test results 0.079AC and 0.079 AC.
- () 10. Air blank completed.
- () 11. When test complete, remove printed record. Attach the record to the completed checklist.
- () 12. Return mode selection switch #9 to off position after all calibration checks are complete or type Q and enter on keyboard.

SIGNATURE *Terry Gallegos*

*WITH OR WITHOUT VAPOR RECIRCULATION AND WITH OR WITHOUT KEYBOARD

DPS Form Exh E-2 (Rev. 05-1)

6076
TSH
jli

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-0055304
09/08/05 16:46

TEST BRAC VALUE TIME
AIR BLANK .000 16:43
*SUBJECT TEST .254 16:46
AIR BLANK .000 16:46

* DEFICIENT SAMPLE - VALUE
PRINTED WAS HIGHEST OBTAINED.

SUBJECT'S NAME
DUI HQ
INSTRUMENT LOCATION

TIME FIRST OBSERVED

OPERATOR
Gallegos 20847

ADDITIONAL INFORMATION AND/OR REMARKS



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.

TUCSON POLICE DEPARTMENT
INTOXILYZER - ALCOHOL ANALYZER
AZ MODEL 5000 SN 66-0055304
09/08/06 16:42

TEST BRAC VALUE TIME
AIR BLANK .000 16:41
INVALID SAMPLE .NDA 16:42
AIR BLANK .000 16:42

BEGIN NEW DEPRIVATION PERIOD

SUBJECT'S NAME
DUI HQ
INSTRUMENT LOCATION

TIME FIRST OBSERVED

OPERATOR
Gallegos 20847

ADDITIONAL INFORMATION AND/OR REMARKS

Sweet breath



INTOXILYZER® INSTRUMENT PRINTER CARD

© 1986 by CMI INC.