



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

B6  
8.18.2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-18-11 Cal 1 0.095 AC Cal 2 0.095 AC

Between the dates of 7.21.11 and 8-18-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


8-18-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 18th, day of August, 20 11.

Susan A. Hitchye

Notary Public

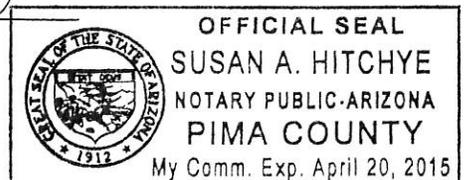

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1067  
CW

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

BH  
8.18.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 8-18-11 TIME 0711

INTOXILYZER SERIAL # 50-000331 LOCATION ODW DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

( ✓ ) 1. Clock time check.

( ✓ ) 2. Date check.

OPERATIONAL TESTS

( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.

( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/18/2011 07:13:18

2007  
CW  
BH

---

Standard Lot#: 654682  
Last Changed By: C. CONLEY #34517

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 07:14:11 |
| Diagnostic Test | Pass   | 07:14:39 |
| Air Blank       | 0.000  | 07:15:06 |
| 0.100 Cal Check | 0.095  | 07:15:25 |
| Air Blank       | 0.000  | 07:15:56 |
| Subject Test    | 0.000  | 07:17:01 |
| Air Blank       | ABT*   | 07:17:08 |
| Air Blank       | 0.000  | 07:17:36 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/18/2011

07:17:55

Standard Lot#: 654682  
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | RFI*   | 07:18:27 |
| Air Blank | IPS**  | 07:18:37 |

\*RFI Detect  
\*\*Improper Sample

34517  
lu  
BH

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/18/2011 07:19:14

4567  
CW  
BY

Standard Lot#: 654682  
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 07:20:00 |
| Diagnostic Test | Pass   | 07:20:27 |
| Air Blank       | 0.000  | 07:20:55 |
| 0.100 Cal Check | 0.095  | 07:21:13 |
| Air Blank       | 0.000  | 07:21:44 |
| Subject Test    | 0.157* | 07:24:52 |
| Air Blank       | 0.000  | 07:25:25 |
| Air Blank       | 0.000  | 07:25:52 |
| Subject Test    | XXX**  | 07:26:22 |
| Air Blank       | 0.000  | 07:26:55 |

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

08/18/2011

07:27:13

-----  
Standard Lot#: 654682

Last Changed By: C. CONLEY #34517

-----  
QAS: C. CONLEY #34517

TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 07:27:41 |
| Cal Check | 0.095  | 07:27:59 |
| Air Blank | 0.000  | 07:28:30 |

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

5087  
CW  
104

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

b.87  
Ch

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
8-18-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 8-18-11 TIME 0711  
Intoxilyzer Serial # 80-000331 Location ODW DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 095 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000331  
08/18/2011 07:11:39

7587  
CW

-----  
Standard Lot#: 654682  
Last Changed By: C. CONLEY #34517  
-----

TS4

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 07:12:07 |
| Cal Check | 0.095  | 07:12:25 |
| Air Blank | 0.000  | 07:12:55 |

*Clare Conley*

\_\_\_\_\_  
Operator's Signature