



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
8-8-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-31-07 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 7-5-07 and 7-31-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

8-3-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 3, day of August, 2007.

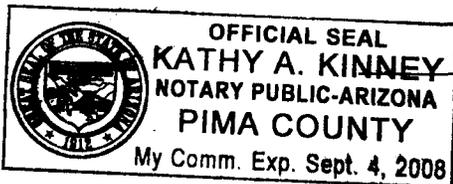

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
8807

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 7-31-07 TIME 10:53

INTOXILYZER SERIAL # 80-000323 LOCATION TPD Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None

SIGNATURE Clare Conley

2 of 8 cu
TSH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
07/31/2007 10:49:31

Standard Lot#: 706602I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:50:23
Diagnostic Test	Pass	10:50:50
Air Blank	0.000	10:51:18
0.100 Cal Check	0.098	10:51:36
Air Blank	0.000	10:52:06
Subject Test	0.000	10:52:23
Air Blank	ABT*	10:52:29
Air Blank	0.000	10:52:58

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 8 cu

Tsh

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000323

07/31/2007

11:00:42

Standard Lot#: 706602I

Last Changed By: T. GALLEGOS

#28847

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:01:36
Diagnostic Test	Pass	11:02:04
Air Blank	RFI*	11:02:18
Air Blank	0.000	11:02:46

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 cu

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
07/31/2007 11:03:23

Standard Lot#: 706602I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	11:03:57
Air Blank	0.000	11:04:26

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 8 cu

Bh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
07/31/2007

11:04:30

Standard Lot#: 706602I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:05:22
Diagnostic Test	Pass	11:05:50
Air Blank	0.000	11:06:17
0.100 Cal Check	0.097	11:06:36
Air Blank	0.000	11:07:05
Subject Test	0.121*	11:10:13
Air Blank	0.000	11:10:46
Air Blank	0.000	11:11:13
Subject Test	XXX**	11:11:32
Air Blank	0.000	11:12:05

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
07/31/2007 11:12:34

6 of 8 cu

BM

Standard Lot#: 706602I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:13:02
Cal Check	0.097	11:13:20
Air Blank	0.000	11:13:50

Clare Conley

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8 cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

sh

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 7-31-07 TIME 10:53
Intoxilyzer Serial # 80-000323 Location TPD Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

8 of 8 u

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
07/31/2007 10:53:24

TS4

Standard Lot#: 706602I
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:53:52
Cal Check	0.098	10:54:10
Air Blank	0.000	10:54:40

Clae Conley

Operator's Signature