



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS
8.3.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/29/10 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 07/01/10 and 07/29/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

08/02/10


ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 2, day of August, 20 10.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

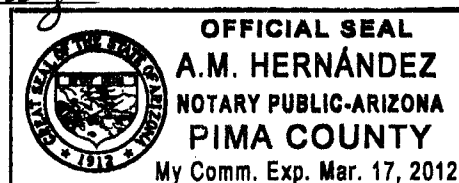

187
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

154
8.3.2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 07/29/10 TIME 14:31
INTOXILYZER SERIAL # 80-000323 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Cheryl Straub

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

287
CS

TS4
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
07/29/2010

14:33:16

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:34:13
Diagnostic Test	Pass	14:34:40
Air Blank	0.000	14:35:08
0.100 Cal Check	0.096	14:35:26
Air Blank	0.000	14:35:56
Subject Test	0.000	14:36:11
Air Blank	ABT*	14:36:18
Air Blank	0.000	14:36:46

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

By

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
07/29/2010 14:36:51

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

OAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:37:24
Air Blank	IPS**	14:37:31

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

487
CS

B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
07/29/2010

14:39:39

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:40:30
Diagnostic Test	Pass	14:40:58
Air Blank	0.000	14:41:25
0.100 Cal Check	0.096	14:41:44
Air Blank	0.000	14:42:13
Subject Test	0.142*	14:45:21
Air Blank	0.000	14:45:54
Air Blank	0.000	14:46:21
Subject Test	XXX**	14:46:41
Air Blank	0.000	14:47:14

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

537
08

TSY

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

07/29/2010

14:47:32

Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:47:59
Cal Check	0.095	14:48:18
Air Blank	0.000	14:48:47



Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 07/29/10 TIME 14:31
Intoxilyzer Serial # 80-000323 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000323

07/29/2010

14:31:03

777
CS
p4


Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:31:30
Cal Check	0.096	14:31:49
Air Blank	0.000	14:32:18

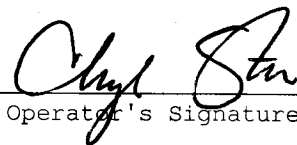
Operator's Signature