



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

CS  
9-27-10

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-23-2010 Cal 1 0.096 AC Cal 2 0.094 AC

Between the dates of 8-26-2010 and 9-23-2010 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
The instrument was returned to service on                     .

Terry Gallegos

Terry S. Gallegos


9-27-2010

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 27, day of September, 2010.

Amy Hernandez

Notary Public

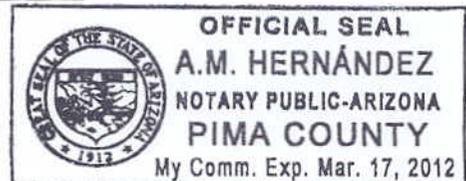

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10f7  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 9-23-2010 TIME 1232

INTOXILYZER SERIAL # 80-000323 LOCATION Crime Lab

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.     
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed   
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 094 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Terry Gallegos

CS

20f7  
B4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
09/23/2010 12:34:42

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:35:35
Diagnostic Test	Pass	12:36:02
Air Blank	0.000	12:36:30
0.100 Cal Check	0.095	12:36:49
Air Blank	0.000	12:37:18
Subject Test	0.000	12:37:36
Air Blank	ABT*	12:37:45
Air Blank	0.000	12:38:14

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

CS

30f7  
TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
09/23/2010 12:38:29

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:39:21
Diagnostic Test	Pass	12:39:49
Air Blank	RFI*	12:39:59
Air Blank	IPS**	12:40:08

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

CS

40F7  
TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000323  
09/23/2010

12:40:24

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:41:16
Diagnostic Test	Pass	12:41:44
Air Blank	0.000	12:42:11
0.100 Cal Check	0.094	12:42:30
Air Blank	0.000	12:42:59
Subject Test	0.168*	12:46:07
Air Blank	0.000	12:46:39
Air Blank	0.000	12:47:07
Subject Test	XXX**	12:47:32
Air Blank	0.000	12:48:06

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
09/23/2010 12:48:33

5077  
TS4

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	12:49:00
Cal Check	0.094	12:49:19
Air Blank	0.000	12:49:48

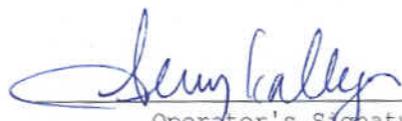
  
-----  
Operator's Signature

EXHIBIT G-2

CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

60f7  
TS4

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 9.23.2010 TIME 1232  
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000323  
09/23/2010 12:32:18

7 of 7  
TS4

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	12:32:45
Cal Check	0.096	12:33:04
Air Blank	0.000	12:33:33

  
-----  
Operator's Signature