



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000323 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

SC
6/10/11

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-2-11 Cal 1 0.095 AC Cal 2 0.094 AC

Between the dates of 5-5-11 and 6-2-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


6-3-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 3RD, day of June, 2011.

Isabel C. Trujillo

Notary Public

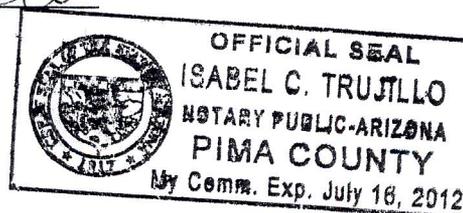

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1067
CM

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
6/10/11

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 6-2-11 TIME 1046
INTOXILYZER SERIAL # 80-000 323 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 094 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None

SIGNATURE Clare Conley

2067
cu

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
06/02/2011

10:47:38

5C

Standard Lot#: 642845
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:48:23
Diagnostic Test	Pass	10:48:51
Air Blank	0.000	10:49:18
0.100 Cal Check	0.094	10:49:37
Air Blank	0.000	10:50:06
Subject Test	0.000	10:50:31
Air Blank	ABT*	10:50:37
Air Blank	0.000	10:51:06

*Sequence Aborted

Not a Successfully Completed Test Sequence

347
CW

SC

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
06/02/2011 10:51:11

Standard Lot#: 642845
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:51:59
Diagnostic Test	Pass	10:52:26
Air Blank	RFI*	10:52:34
Air Blank	IPS**	10:52:41

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000323
06/02/2011

10:52:47

4567
CW

5C

Standard Lot#: 642845
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:53:35
Diagnostic Test	Pass	10:54:03
Air Blank	0.000	10:54:30
0.100 Cal Check	0.095	10:54:49
Air Blank	0.000	10:55:18
Subject Test	0.083*	10:58:26
Air Blank	0.000	10:58:58
Air Blank	0.000	10:59:25
Subject Test	XXX**	10:59:45
Air Blank	0.000	11:00:18

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
06/02/2011 11:00:43

587
CN

SC

Standard Lot#: 642845
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:01:10
Cal Check	0.094	11:01:29
Air Blank	0.000	11:01:58

Clare Conley

Operator's Signature

EXHIBIT G-2

6567
CW

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

SC
6/10/11

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 6-2-11 TIME 1046
Intoxilyzer Serial # 80-000323 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000323
06/02/2011 10:46:07

7567
en

SC

Standard Lot#: 642845
Last Changed By: T. GALLEGOS #28847

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:46:35
Cal Check	0.095	10:46:54
Air Blank	0.000	10:47:23

Clare Conley

Operator's Signature