



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

By  
8-29-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-29-08 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 7-31-08 and 8-29-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

8-29-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 29<sup>th</sup> day of August, 2008.

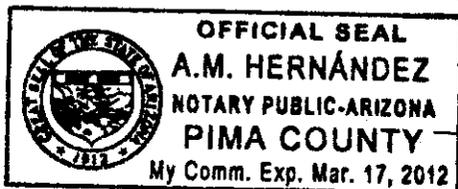

A.M. Hernandez

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

1067  
On  
154  
8-29-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 8-29-08 TIME 0650  
INTOXILYZER SERIAL # 80-000331 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/29/2008

06:52:21

2 of 7  
C  
B4

Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:53:11
Diagnostic Test	Pass	06:53:39
Air Blank	0.000	06:54:06
0.100 Cal Check	0.098	06:54:25
Air Blank	0.000	06:54:54
Subject Test	0.000	06:55:29
Air Blank	ABT*	06:55:35
Air Blank	0.000	06:56:03

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/29/2008

06:56:19

387  
cu  
PM

Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:56:50
Air Blank	IPS**	06:56:59

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000331  
08/29/2008

06:57:18

487  
er  
BY

Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:58:08
Diagnostic Test	Pass	06:58:36
Air Blank	0.000	06:59:03
0.100 Cal Check	0.098	06:59:22
Air Blank	0.000	06:59:52
Subject Test	0.204*	07:03:01
Air Blank	0.000	07:03:33
Air Blank	0.000	07:04:00
Subject Test	XXX**	07:04:22
Air Blank	0.000	07:04:55

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000331  
08/29/2008 07:06:00

-----  
Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:06:27
Cal Check	0.098	07:06:46
Air Blank	0.000	07:07:16

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

5567  
C  
- 1/4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

6587  
cu  
BY  
8-29-08

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 8-29-08 TIME 0650  
Intoxilyzer Serial # 80-000331 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000331  
08/29/2008 06:50:37

-----  
Standard Lot#: 719801I  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:51:04
Cal Check	0.098	06:51:23
Air Blank	0.000	06:51:52

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

7887  
cu  
BY



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

TSY  
8.29.08

Date: 8-29-08  
QAS: Clare Conley  
Intoxilyzer Serial Number: 80-000331  
New Lot Number: 811201 I  
New Expiration Date: 31 Mar 2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0.101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000331  
08/29/2008 07:19:22

25/2

cu

154

-----  
Standard Lot#: 811201I  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:19:49
Cal Check	0.101	07:20:08
Air Blank	0.000	07:20:37

*Clare Conley*

\_\_\_\_\_  
Operator's Signature