



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
6-10-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000331 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 06/04/09 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 05/07/09 and 06/04/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.

Cheryl K. Straub

Cheryl K. Straub


06/09/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 09th, day of June, 2009.

A.M. Hernández

Notary Public

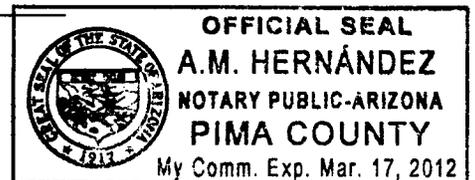

*Includes documents specified in A.A.C. R13-10-104.A.5

189
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

By
6-10-09

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 06/04/09 TIME 08:46
INTOXILYZER SERIAL # 80-000331 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS Alcohol free subject test repeated due
to power interruption; intox and printer
were undrugged during operation.

SIGNATURE Cheryl Straub

279
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
06/04/2009

08:49:04

ADY

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:49:48
Air Blank	IPS**	08:50:07

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

379
C8

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
06/04/2009 08:50:37

TS4

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:51:37
Diagnostic Test	Pass	08:52:05
Air Blank	0.000	08:52:32
0.100 Cal Check	0.100	08:52:51
Air Blank	0.000	08:53:21
Subject Test	0.000	08:53:39
Air Blank	ABT*	08:53:46
Air Blank	0.000	08:54:15

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

489
08

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
06/04/2009

09:05:55

104

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:07:11
Diagnostic Test	Pass	09:07:39
Air Blank	0.000	09:08:06
0.100 Cal Check	0.100	09:08:25
Air Blank	0.000	09:08:55
Subject Test	0.000	09:09:14
Air Blank	ABT*	09:09:21
Air Blank	0.000	09:09:49

*Sequence Aborted

TEST REPEATED DUE TO POWER INTERRUPTION
08-50772

Not a Successfully Completed Test Sequence

589
CS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
06/04/2009

09:12:39

D4

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:13:41
Diagnostic Test	Pass	09:14:08
Air Blank	0.000	09:14:36
0.100 Cal Check	0.100	09:14:54
Air Blank	0.000	09:15:25
Subject Test	XXX*	09:16:53
Air Blank	0.000	09:17:26

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

689
CS

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000331
06/04/2009 10:45:38

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:46:48
Diagnostic Test	Pass	10:47:15
Air Blank	0.000	10:47:43
0.100 Cal Check	0.100	10:48:02
Air Blank	0.000	10:48:32
Subject Test	0.051*	10:51:40
Air Blank	0.000	10:52:12
Air Blank	ABT**	10:52:22
Air Blank	0.000	10:52:50

*Deficient Sample
**Sequence Aborted

Not a Successfully Completed Test Sequence

789
CS


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000331
06/04/2009 09:18:47

TSM

Standard Lot#: 811201I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:19:15
Cal Check	0.099	09:19:33
Air Blank	0.000	09:20:03

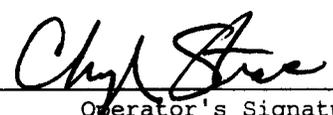

Operator's Signature

EXHIBIT G-2

889
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TJ4
6-10-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 06/04/09 TIME 08:46
Intoxilyzer Serial # 80-000331 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

999
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000331

06/04/2009

08:46:13

154


Standard Lot#: 811201I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:46:40
Cal Check	0.100	08:46:59
Air Blank	0.000	08:47:29

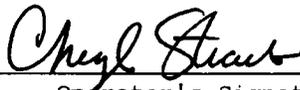

Operator's Signature