



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

2587

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4.13.2007 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 3.15.2007 and 4.13.2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.


Terry S. Gallegos

Terry S. Gallegos

4.24.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24, day of April, 2007.

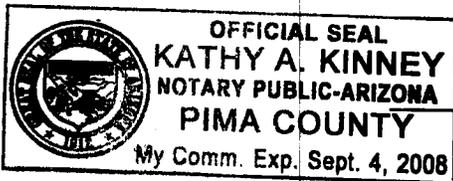

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10 of 7
D/MP

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 4.13.2007 TIME 1811
INTOXILYZER SERIAL # 80-000332 LOCATION ODE

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning. ✓ ✓
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
04/13/2007

18:13:19

2 of 7
B6mm

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 18:14:11 |
| Diagnostic Test | Pass | 18:14:38 |
| Air Blank | 0.000 | 18:15:06 |
| 0.100 Cal Check | 0.100 | 18:15:24 |
| Air Blank | 0.000 | 18:15:54 |
| Subject Test | 0.000 | 18:16:14 |
| Air Blank | 0.000 | 18:16:42 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 18:17:17 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
04/13/2007

18:17:22

3077
Bhmm

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 18:18:12 |
| Diagnostic Test | Pass | 18:18:39 |
| Air Blank | RFI* | 18:18:51 |
| Air Blank | IPS** | 18:19:12 |

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
04/13/2007

18:19:23

4077
TSGM

Standard Lot#: 630601I
Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 18:20:17 |
| Diagnostic Test | Pass | 18:20:44 |
| Air Blank | 0.000 | 18:21:11 |
| 0.100 Cal Check | 0.100 | 18:21:30 |
| Air Blank | 0.000 | 18:22:00 |
| Subject Test | 0.021* | 18:25:08 |
| Air Blank | 0.000 | 18:25:39 |
| Air Blank | 0.000 | 18:26:07 |
| Subject Test | XXX** | 18:26:28 |
| Air Blank | 0.000 | 18:27:01 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

04/13/2007

18:27:26

*SofS
D/yr*

Standard Lot#: 630601I

East Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test g/210L Time

Air Blank 0.000 18:27:53
Cal Check 0.100 18:28:12
Air Blank 0.000 18:28:42



Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6077
TJG

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallenas 200547 AGENCY Tucson Police Department
DATE 4.13.2007 TIME 1911
Intoxilyzer Serial # 80-000332 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallenas

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

04/13/2007

18:11:56

7 of 7
TS, MD

Standard Lot#: 6306011

Last Changed By: J. COOMBS #46160

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test g/210L Time

Air Blank 0.000 18:12:23
Cal Check 0.100 18:12:42
Air Blank 0.000 18:13:11



Operator's Signature