



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
8-1-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-25-08 Cal 1 0.098 AC Cal 2 0.096 AC

Between the dates of 7-3-08 and 7-25-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A
The instrument was returned to service on N/A


Clare Conley

Clare L. Conley

7-25-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 25th day of July, 2008.

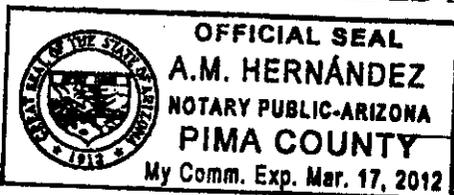

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1569
Cu
754
8-1-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 7-25-08 TIME 0734

INTOXILYZER SERIAL # 80-000332 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
07/25/2008

07:36:36

269
Cm
P4

Standard Lot#: 735401G
Last Changed By: T. O HARA #38876

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:37:24
Diagnostic Test	Pass	07:37:52
Air Blank	0.000	07:38:19
0.100 Cal Check	0.097	07:38:38
Air Blank	0.000	07:39:08
Subject Test	0.000	07:39:23
Air Blank	ABT*	07:39:29
Air Blank	0.000	07:39:58

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
07/25/2008

07:43:26

3569
an
by

Standard Lot#: 735401G
Last Changed By: T. O HARA #38876

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:43:57
Air Blank	0.000	07:44:26

*RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
07/25/2008

07:44:35

469
Cu
PH

Standard Lot#: 735401G
Last Changed By: T. O HARA #38876

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:45:34
Diagnostic Test	Pass	07:46:01
Air Blank	0.000	07:46:28
0.100 Cal Check	0.097	07:46:47
Air Blank	0.000	07:47:17
Subject Test	0.121*	07:50:25
Air Blank	0.000	07:50:58
Air Blank	0.000	07:51:25
Subject Test	XXX**	07:51:54
Air Blank	0.000	07:52:28

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
07/25/2008 07:52:46

5 of 9
cu
104

Standard Lot#: 735401G
Last Changed By: T. O HARA #38876

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:53:14
Cal Check	0.096	07:53:32
Air Blank	0.000	07:54:02

Clare Conley

Operator's Signature

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
07/25/2008 09:17:20

6069
Cu
754

Standard Lot#: 735401G
Last Changed By: T. O HARA #38876

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	09:18:04
Air Blank	0.000	09:18:32

*Improper Sample

Not a Successfully Completed Test Sequence

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7569
cu
756
8/1/08

QA Specialist Clare Conley AGENCY Tucson Police Department

DATE 7-25-08 TIME 0734

Intoxilyzer Serial # 80-000332 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 150 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

07/25/2008

07:34:03

Standard Lot#: 735401G

Last Changed By: T. O HARA #38876

QAS: C. CONLEY #34517

TUCSON PD

Test g/210L Time

Air Blank 0.000 07:34:30
Reference ABT* 07:34:38
Air Blank 0.000 07:35:06

*Sequence Aborted

Clare Conley

Operator's Signature

809

cu


104

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
07/25/2008 07:35:12

Standard Lot#: 735401G
Last Changed By: T. O HARA #38876

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:35:39
Cal Check	0.098	07:35:58
Air Blank	0.000	07:36:28

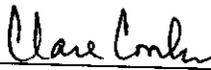

Operator's Signature

9069
cu

104