



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
9-25-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000 332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-19-08 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 8-25-08 and 9-19-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley
Clare L. Conley

9-25-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 25th day of September, 2008.

A. Hernandez
Notary Public

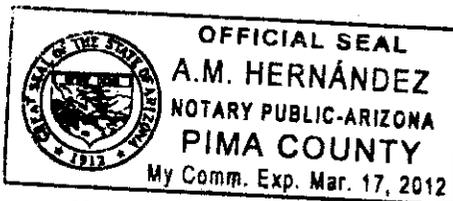

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1067
en
754
9-25-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 9-19-08 TIME 0739
INTOXILYZER SERIAL # 80-000332 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 005 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
09/19/2008

07:54:23

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 07:55:14 |
| Diagnostic Test | Pass | 07:55:41 |
| Air Blank | 0.000 | 07:56:08 |
| 0.100 Cal Check | 0.097 | 07:56:27 |
| Air Blank | 0.000 | 07:56:57 |
| Subject Test | 0.000 | 07:57:13 |
| Air Blank | ABT* | 07:57:20 |
| Air Blank | 0.000 | 07:57:48 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

287
cu
B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
09/19/2008 07:59:15

347
an
ph

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 07:59:42 |
| Air Blank | IPS** | 07:59:49 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
09/19/2008

08:00:34

Standard Lot#: 735401G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:01:23 |
| Diagnostic Test | Pass | 08:01:50 |
| Air Blank | 0.000 | 08:02:18 |
| 0.100 Cal Check | 0.096 | 08:02:37 |
| Air Blank | 0.000 | 08:03:06 |
| Subject Test | 0.204* | 08:06:15 |
| Air Blank | 0.000 | 08:06:48 |
| Air Blank | 0.000 | 08:07:15 |
| Subject Test | XXX** | 08:07:34 |
| Air Blank | 0.000 | 08:08:07 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

487
en
B4

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

09/19/2008

08:08:26

Standard Lot#: 735401G

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 08:08:53 |
| Cal Check | 0.095 | 08:09:12 |
| Air Blank | 0.000 | 08:09:42 |

Clare Conley

Operator's Signature

5067
ew
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6567
TBH
9.25.08
LW

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 9-19-08 TIME 0739
Intoxilyzer Serial # 80-100332 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Qfi

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

09/19/2008

07:39:08

Standard Lot#: 735401G

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 07:39:36 |
| Cal Check | 0.097 | 07:39:54 |
| Air Blank | 0.000 | 07:40:24 |

Claire Conley

Operator's Signature

7887
Qu
124

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 2 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TSY
10-2-2008

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/29/08 TIME 10:21
Intoxilyzer Serial # 80-000332 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 2 sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
09/29/2008 10:21:15

Standard Lot#: 7354016
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 10:21:42 |
| Cal Check | 0.097 | 10:22:01 |
| Air Blank | 0.000 | 10:22:31 |

Sarah Cione
Operator's Signature