

1 of 85C

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

134
2-18-09

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 1/14/09 TIME 15:17

INTOXILYZER SERIAL # 80-000332 LOCATION Vehicle #617 - Team 5

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 8 sc

TS4

INTOXILYZER 6000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
01/14/2009 15:18:37

Standard Lot#: 808101G
Last Changed By: C. ANDREACOLA #27359

DAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:19:23
Diagnostic Test	Pass	15:19:52
Air Blank	0.000	15:20:19
0.100 Cal Check	0.098	15:20:38
Air Blank	0.000	15:21:08
Subject Test	0.000	15:21:25
Air Blank	0.000	15:21:53
Five Minute Wait	ABT*	
Air Blank	0.000	15:22:28

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 8 sc

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
01/14/2009 15:22:45

Standard Lot#: 808101G
Last Changed By: C. ANDREACOLA #27359

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	15:23:25
Air Blank	IPS**	15:23:36

*RFI Detect

**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 8 sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
01/14/2009 15:24:54

Standard Lot#: 808101G
Last Changed By: C. ANDREACOLA #27359

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:25:57
Diagnostic Test	Pass	15:26:26
Air Blank	0.000	15:26:53
0.100 Cal Check	0.099	15:27:12
Air Blank	0.000	15:27:41
Subject Test	0.109*	15:30:50
Air Blank	0.000	15:31:22
Air Blank	0.000	15:31:50
Subject Test	XXX**	15:32:10
Air Blank	0.000	15:32:44

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 8 SC

TJH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
01/14/2009 15:33:21

Standard Lot#: 808101G
Last Changed By: C. ANDREACOLA #27359

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:33:48
Cal Check	0.099	15:34:07
Air Blank	RFI*	15:34:17
Air Blank	RFI*	15:34:25

*RFI Detect

Sarah Cione
Operator's Signature

6 of 85c

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
01/14/2009 15:35:30

Standard Lot#: 808101G
Last Changed By: C. ANDREACOLA #27359

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:35:58
Cal Check	0.098	15:36:16
Air Blank	0.000	15:36:46

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

134
2-18-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department

DATE 1/14/09 TIME 15:17

Intoxilyzer Serial # 80-000332 Location Vehicle # 617 - Team 5

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

80785C

154

INTOXILYZER 800G
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
01/14/2009 15:17:11

Standard Lot#: 808101G
Last Changed By: C. ANDREACOLA #27359

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:17:38
Cal Check	0.098	15:17:57
Air Blank	0.000	15:18:27

Sarah Cione
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

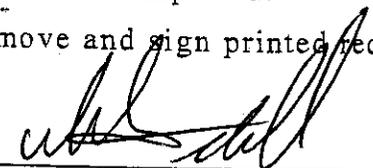
1 of 2 sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

TS4
2.18.09

Date: 12/19/2008
QAS: William D. Horanick
Intoxilyzer Serial Number: 80-000332
New Lot Number: 808101G
New Expiration Date: 3/31/2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 0189 AC.
TS4 2.18.09
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE  

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
12/19/2008

16:42:39

2 of 25c

Standard Lot#: 808101G
Last Changed By: J. TRAINING #99999

By
2-18-09

QAS: J. TRAINING #99999
AZ DPS

Test	g/210L	Time
Air Blank	0.000	16:43:06
Cal Check	0.099	16:43:25
Air Blank	0.000	16:43:54


Operator's Signature 