



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*Handwritten:* 7-1-09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-25-09 Cal 1 0.098 AC Cal 2 0.095 AC

Between the dates of 6-5-09 and 6-25-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


6-30-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 30<sup>th</sup>, day of June, 2009.

A.M. Hernandez

Notary Public

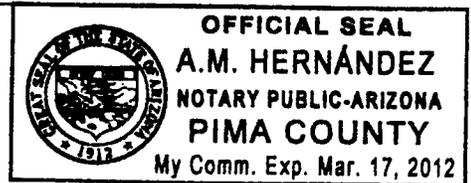

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
cu  
BY  
7.1.09

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 6-25-09 TIME 0930

INTOXILYZER SERIAL # 80-000332 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000332  
06/25/2009 09:31:54

287  
en

Standard Lot#: 811201I  
Last Changed By: S. CIONE #50616

BH

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:32:48
Diagnostic Test	Pass	09:33:15
Air Blank	0.000	09:33:43
0.100 Cal Check	0.098	09:34:02
Air Blank	0.000	09:34:31
Subject Test	0.000	09:34:49
Air Blank	ABT*	09:34:55
Air Blank	0.000	09:35:24

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.*

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000332  
06/25/2009 09:35:33

367  
Cu  
154

---

Standard Lot#: 811201I  
Last Changed By: S. CIONE #50616

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	09:36:06
Air Blank	IPS**	09:36:14

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.*

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000332  
06/25/2009

09:37:53

4567  
cu

Standard Lot#: 811201I  
Last Changed By: S. CIONE #50616

104

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:38:39
Diagnostic Test	Pass	09:39:06
Air Blank	0.000	09:39:34
0.100 Cal Check	0.097	09:39:52
Air Blank	0.000	09:40:22
Subject Test	0.439*	09:43:30
Air Blank	0.000	09:44:05
Air Blank	0.000	09:44:32
Subject Test	XXX**	09:45:15
Air Blank	0.000	09:45:50

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

06/25/2009

09:46:08

-----  
Standard Lot#: 8112011

Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:46:36
Cal Check	0.095	09:46:54
Air Blank	0.000	09:47:24

587  
Cu  
B4

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

657  
CM  
TBY  
7-1-09

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 6-25-09 TIME 0930  
Intoxilyzer Serial # 80-000332 Location ODW DU HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000332

06/25/2009

09:30:27

-----  
Standard Lot#: 811201I

Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:30:54
Cal Check	0.098	09:31:13
Air Blank	0.000	09:31:43

*C. Conley*

\_\_\_\_\_  
Operator's Signature

7 of 7  
Cn

ASg