



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS  
55-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/29/10 Cal 1 0.097 AC Cal 2 0.099 AC

Between the dates of 4/1/10 and 4/29/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/30/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30<sup>th</sup>, day of April, 2010.

A.M. Hernández

Notary Public

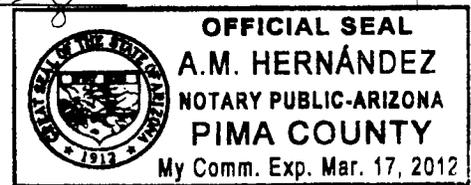

\*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 8 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY  
6.5.2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/29/10 TIME 11:01

INTOXILYZER SERIAL # 80-000332 LOCATION Vehicle #945

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

2 of 8 sc

104

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000332  
04/29/2010 11:02:31

-----  
Standard Lot#: 642845  
Last Changed By: C. CONLEY #34517  
-----

OAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:03:23
Diagnostic Test	Pass	11:03:52
Air Blank	0.000	11:04:19
0.100 Cal Check	0.097	11:04:38
Air Blank	0.000	11:05:08
Subject Test	0.000	11:05:22
Air Blank	0.000	11:05:50
Five Minute Wait	ABT*	
Air Blank	0.000	11:06:25

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 8 sc

134

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000332  
04/29/2010 11:06:31

Standard Lot#: 642845  
Last Changed By: C. CONLEY #34517

DAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:06:59
Air Blank	IPS**	11:07:09

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 8 SC

134

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000332  
04/29/2010 11:07:16

Standard Lot#: 642845  
Last Changed By: C. CONLEY #34517

OAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:08:02
Diagnostic Test	Fail*	11:08:33
Air Blank	0.000	11:09:01

\* Temperature Regulation Test

Not a Successfully  
Completed  
Test Sequence

5 of 8 sc

TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000332  
04/29/2010 11:09:36

Standard Lot#: 642845  
Last Changed By: C. CONLEY #34517

OAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:10:43
Diagnostic Test	Pass	11:11:12
Air Blank	0.000	11:11:39
0.100 Cal Check	0.098	11:11:58
Air Blank	0.000	11:12:28
Subject Test	0.070*	11:15:36
Air Blank	0.000	11:16:08
Air Blank	0.000	11:16:35
Subject Test	XXX**	11:16:54
Air Blank	0.000	11:17:27

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

6 of 8 sc

T04

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000332  
04/29/2010 11:17:54

-----  
Standard Lot#: 642845  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:18:21
Cal Check	0.099	11:18:40
Air Blank	0.000	11:19:10

*Sarah Cione*  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

7 of 8 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

B4  
5.5.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department

DATE 4/29/10 TIME 11:01

Intoxilyzer Serial # 80-000332 Location Vehicle #945

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

P of 8 sc

134

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000332  
04/29/2010 11:01:08

-----  
Standard Lot#: 642845  
Last Changed By: C. CONLEY #34517  
-----

OAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:01:35
Cal Check	0.097	11:01:54
Air Blank	0.000	11:02:24

*Sarah Cione*  
-----  
Operator's Signature