



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
1-18-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/5/12 Cal 1 0.095 AC Cal 2 0.096 AC

Between the dates of 12/8/11 and 1/5/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


1/11/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 11, day of January, 2012.

A.M. Hernandez

Notary Public

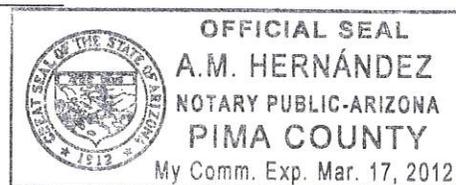

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7sc

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TS4
1.18.2012

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 1/5/12 TIME 1354
INTOXILYZER SERIAL # 80-000332 LOCATION ODW-Tm 2

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7_{sc}

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
01/05/2012

13:55:56

134

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:56:38
Diagnostic Test	Pass	13:57:06
Air Blank	0.000	13:57:33
0.100 Cal Check	0.096	13:57:52
Air Blank	0.000	13:58:22
Subject Test	0.000	13:58:38
Air Blank	0.000	13:59:05
Five Minute Wait	ABT*	
Air Blank	0.000	13:59:40

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7 sc

B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
01/05/2012 13:59:46

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:00:38
Air Blank	IPS**	14:00:46

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000332
01/05/2012

14:00:53

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:01:39
Diagnostic Test	Pass	14:02:06
Air Blank	0.000	14:02:33
0.100 Cal Check	0.097	14:02:52
Air Blank	0.000	14:03:22
Subject Test	0.044*	14:06:30
Air Blank	0.000	14:07:03
Air Blank	0.000	14:07:30
Subject Test	XXX**	14:07:56
Air Blank	0.000	14:08:31

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
01/05/2012 14:08:48

TS4

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:09:16
Cal Check	0.096	14:09:34
Air Blank	0.000	14:10:04



Operator's Signature

6 of 7sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
1.18.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 1/5/12 TIME 1354
Intoxilyzer Serial # 80-000332 Location ODW-Tm 2

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000332
01/05/2012 13:54:23

7 of 7sc

134

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:54:50
Cal Check	0.095	13:55:09
Air Blank	0.000	13:55:39

Sarah Cione

Operator's Signature