



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
7/31/12

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7.26.2012 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 6.28.2012 and 7.26.2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on [blank]
The instrument was returned to service on [blank]

[Handwritten signature of Terry S. Gallegos]

Terry S. Gallegos


7.31.2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 31, day of July, 2012.

[Handwritten signature of Ann L. Beecroft]

Notary Public

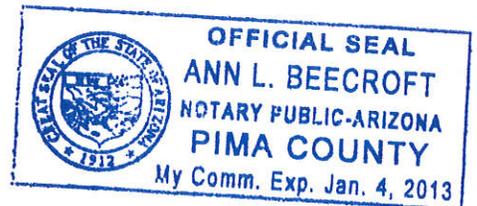

\*Includes documents specified in A.A.C. R13-10-104.A.5

10 of 7  
154

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

SC  
7/31/12

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 7-26-2012 TIME 1153

INTOXILYZER SERIAL # 80-000332 LOCATION ODS

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7  
154

SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/26/2012 11:56:42

---

Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:57:35
Diagnostic Test	Pass	11:58:10
Air Blank	0.000	11:58:38
0.100 Cal Check	0.097	11:58:59
Air Blank	0.000	11:59:28
Subject Test	0.000	11:59:45
Air Blank	0.000	12:00:12
Wait	ABT*	
Air Blank	0.000	12:00:48

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7  
TS4

SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/26/2012 12:00:54

---

Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:01:46
Diagnostic Test	Pass	12:02:21
Air Blank	0.000	12:02:48
0.100 Cal Check	0.097	12:03:09
Air Blank	RFI*	12:03:16
Air Blank	IPS**	12:03:24

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/26/2012

12:03:34

4 of 7  
TS4  
SC

---

Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:04:27
Diagnostic Test	Pass	12:05:02
Air Blank	0.000	12:05:30
0.100 Cal Check	0.097	12:05:50
Air Blank	0.000	12:06:20
Subject Test	0.128*	12:09:28
Air Blank	0.000	12:10:00
Air Blank	0.000	12:10:27
Subject Test	XXX**	12:10:41
Air Blank	0.000	12:11:16

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/26/2012 12:11:32

507  
TS4

-----  
Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:12:03
0.100 Cal Check	0.097	12:12:24
Air Blank	0.000	12:12:53

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
TJh

SC  
7/31/12

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 7.26.2012 TIME 1153  
Intoxilyzer Serial # 80-000332 Location ODS

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-000332  
Core Version: 8105.48  
07/26/2012 11:53:22

7 of 7  
TS4

-----  
Standard Lot#: 04911100A1  
Last Changed By: S. CIONE #50616  
-----

SC

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:54:28
0.100 Cal Check	0.098	11:54:49
Air Blank	0.000	11:55:19

Preventative Maintenance  
Performed -- Timer Reset