



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000332 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

SC
8/30/12

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-23-2012 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 7-26-2012 and 8-23-2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry Gallegos

Terry S. Gallegos


8-28-2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 28th, day of August, 2012.

Elma Escarcega

Notary Public

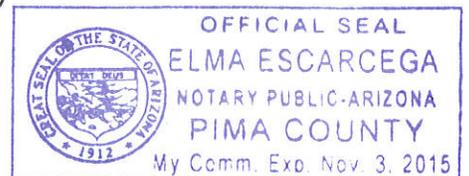

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10f9
TSh

SC
8/30/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 8.23.2012 TIME 1703

INTOXILYZER SERIAL # 80-000332 LOCATION ODS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS 1.75 hr break for court

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
08/23/2012

17:03:17

20f9
TS4

SC

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:04:12
Diagnostic Test	Pass	17:04:47
Air Blank	0.000	17:05:14
0.100 Cal Check	0.099	17:05:35
Air Blank	0.000	17:06:05
Subject Test	0.000	17:06:24
Air Blank	0.000	17:06:52
Wait	ABT*	
Air Blank	0.000	17:07:27

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
08/23/2012

17:07:33

30f 9
TS4

5C

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:08:36
Diagnostic Test	Pass	17:09:11
Air Blank	RFI*	17:09:29
Air Blank	IPS**	17:09:40

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
08/23/2012

17:09:44

4079
B4

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:10:36
Diagnostic Test	Pass	17:11:11
Air Blank	0.000	17:11:38
0.100 Cal Check	0.099	17:11:59
Air Blank	0.000	17:12:29
Subject Test	0.154*	17:15:36
Air Blank	0.000	17:16:08
Air Blank	0.000	17:16:36
Subject Test	XXX**	17:16:55
Air Blank	0.000	17:17:29

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
08/23/2012 17:17:59

Sof 9
B4

SC

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:18:30
0.100 Cal Check	0.099	17:18:51
Air Blank	0.000	17:19:21

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

bof 9
B4

SC
8/30/12

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 8-23-2012 TIME 1513
Intoxilyzer Serial # 80-000332 Location ODS

- () 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- () 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- () 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- () 4. Air blank completed.
- () 5. Calibration check completed. Test results 0. 099 AC.
- () 6. Air blank completed.
- () 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
08/23/2012 15:13:40

7089
B4
SC

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:14:22
0.100 Cal Check	0.099	15:14:43
Air Blank	0.000	15:15:12

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

80F9
TS4

SC
8/30/12

Date: 8-23-2012
 QAS: Terry Gallegos 28847
 Intoxilyzer Serial Number: 80-000332
 New Lot Number: 13712100A2
 New Expiration Date: 7-1-2014

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
 Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
 Select "D" for dry control test. After selection is made press
 ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0.099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and ~~sign printed record~~^{TS4}. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000332
Core Version: 8105.48
08/23/2012 17:21:15

9059
B4
SC

Standard Lot#: 13712100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:21:47
0.100 Cal Check	0.099	17:22:08
Air Blank	0.000	17:22:38

Preventative Maintenance
Performed -- Timer Reset