



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

M  
8-21-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7.19.2007 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 6.21.2007 and 7.19.2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

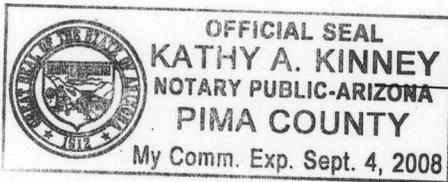
Terry S. Gallegos

Terry S. Gallegos

8.15.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 15, day of August, 2007.



Kathy A. Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

*DU 8/21/07*  
*1 of 9*  
*TJ4*

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 7.19.2007 TIME 1842

INTOXILYZER SERIAL # 80-000334 LOCATION DUI HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed *✓✓✓*
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

*Terry Gallegos*

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
07/19/2007 19:00:27

*mu*  
*2089*  
*184*

---

Standard Lot#: 706602I  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	19:01:21
Diagnostic Test	Pass	19:01:48
Air Blank	0.000	19:02:16
0.100 Cal Check	0.099	19:02:34
Air Blank	0.000	19:03:05
Subject Test	0.000	19:03:41
Air Blank	0.000	19:04:09
Five Minute Wait	ABT*	
Air Blank	0.000	19:04:46

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
07/19/2007

19:05:02

cu

30f 9  
B4

---

Standard Lot#: 706602I  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	19:05:55
Diagnostic Test	Pass	19:06:23
Air Blank	RFI*	19:06:41
Air Blank	IPS**	19:06:51

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
07/19/2007

CM

19:07:01

4 of 9  
BY

---

Standard Lot#: 706602I  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	19:07:53
Diagnostic Test	Pass	19:08:20
Air Blank	0.000	19:08:47
0.100 Cal Check	0.098	19:09:06
Air Blank	0.000	19:09:37
Subject Test	0.104*	19:12:45
Air Blank	0.000	19:13:18
Air Blank	0.000	19:13:45
Subject Test	XXX**	19:14:32
Air Blank	0.000	19:15:07

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

07/19/2007

19:15:55

cm

Sof 9  
Dy

-----  
Standard Lot#: 706602I

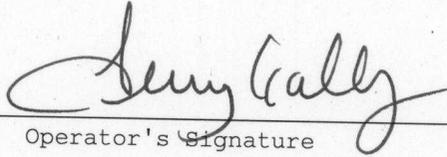
Last Changed By: T. GALLEGOS #28847

-----  
QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

-----  
-----

Test	g/210L	Time
Air Blank	0.000	19:16:22
Cal Check	0.096	19:16:41
Air Blank	0.000	19:17:12



Operator's Signature

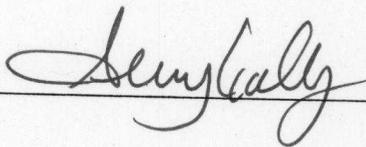
THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

CH 8.21.07  
6 of 9  
TS4

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 7.19.2007 TIME 1842  
Intoxilyzer Serial # 80-000334 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

07/19/2007

18:42:05

-----  
Standard Lot#: 627101G

Last Changed By: T. GALLEGOS #28847

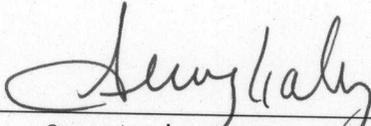
-----  
QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

-----  
-----

Test	g/210L	Time
Air Blank	0.000	18:42:33
Cal Check	0.097	18:42:52
Air Blank	0.000	18:43:22

DU  
7089  
TS4

  
-----  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

CH 8.21-07  
8079

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 7-19-2007  
QAS: Terry Gallegos 20847  
Intoxilyzer Serial Number: 80-000334  
New Lot Number: 706602 I  
New Expiration Date: 3-11-2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

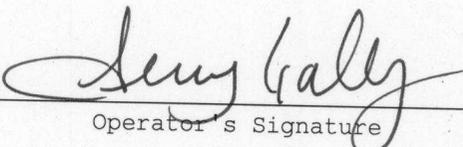
INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
07/19/2007 18:58:10

W  
9 of 9  
TS4

-----  
Standard Lot#: 706602I  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:58:37
Cal Check	0.098	18:58:56
Air Blank	0.000	18:59:26

  
-----  
Operator's Signature