



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

709  
8-29-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-28-08 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 7-31-08 and 8-28-08 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on N/A  
The instrument was returned to service on N/A

Clare Conley

Clare L. Conley

8-29-08

ADPS QAS Permit # 33314

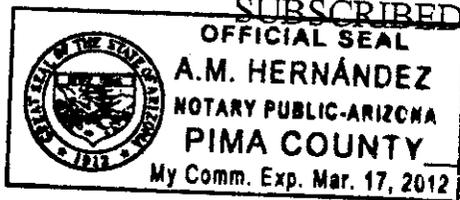

SUBSCRIBED AND SWORN to before me this 29th day of August, 2008.

A.M. Hernandez

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 Cu  
BY 8-29-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 8-28-08 TIME 1125  
INTOXILYZER SERIAL # 80-000334 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
08/28/2008

11:26:36

2 of 7  
CW  
PM

Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:27:24
Diagnostic Test	Pass	11:27:52
Air Blank	0.000	11:28:19
0.100 Cal Check	0.098	11:28:38
Air Blank	0.000	11:29:08
Subject Test	0.000	11:29:24
Air Blank	ABT*	11:29:30
Air Blank	0.000	11:29:59

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
08/28/2008

11:30:30

3567  
eu  
124

---

Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	11:31:00
Air Blank	IPS**	11:31:10

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.*

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
08/28/2008 11:31:53

4507  
cu  
BY

---

Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	11:32:39
Diagnostic Test	Pass	11:33:06
Air Blank	0.000	11:33:33
0.100 Cal Check	0.099	11:33:52
Air Blank	0.000	11:34:22
Subject Test	0.116*	11:37:30
Air Blank	0.000	11:38:02
Air Blank	0.000	11:38:29
Subject Test	XXX**	11:38:44
Air Blank	0.000	11:39:17

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

08/28/2008

11:39:33

-----  
Standard Lot#: 732503I

Last Changed By: S. CIONE #50616

-----  
QAS: C. CONLEY #34517

TUCSON PD

-----  
Test                      g/210L                      Time  
-----  
Air Blank                      0.000                      11:40:01  
Cal Check                      0.099                      11:40:19  
Air Blank                      0.000                      11:40:49

*C. Conley*

\_\_\_\_\_  
Operator's Signature

507  
cu  
py

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

6587  
Cu  
P4 8.29.08

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 8-28-08 TIME 1125  
Intoxilyzer Serial # 80-000334 Location ODW DW HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
08/28/2008 11:25:11

-----  
Standard Lot#: 732503I  
Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	11:25:38
Cal Check	0.099	11:25:57
Air Blank	0.000	11:26:26

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

7507  
ew  
154