



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
12-24-08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-18-08 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 11-20-08 and 12-18-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
 The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


12-23-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 23, day of December, 2008.

A.M. Hernandez

Notary Public

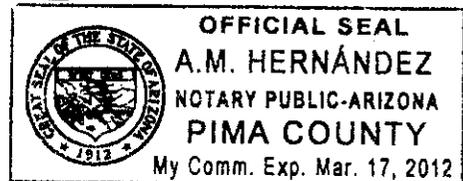

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8
on
TS4 12-24-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 12-18-08 TIME 0938

INTOXILYZER SERIAL # 80-000334 LOCATION DUI HQ ODW

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
12/18/2008 09:39:40

2 of 8

CU

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

TS4

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:40:26
Diagnostic Test	Pass	09:40:53
Air Blank	0.000	09:41:21
0.100 Cal Check	0.100	09:41:39
Air Blank	0.000	09:42:09
Subject Test	0.028	09:42:27
Air Blank	ABT*	09:42:47
Air Blank	0.000	09:43:16

*Sequence Aborted

Repeated alcohol-free subject test.

CU
34517

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
12/18/2008

09:44:02

348
m

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

TS4

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:44:57
Diagnostic Test	Pass	09:45:24
Air Blank	0.000	09:45:52
0.100 Cal Check	0.100	09:46:10
Air Blank	0.000	09:46:40
Subject Test	0.000	09:46:55
Air Blank	ABT*	09:47:02
Air Blank	0.000	09:47:30

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
12/18/2008 09:47:46

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:48:16
Air Blank	IPS**	09:48:24

*RFI Detect
**Improper Sample

4/28
CW

TS4

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
12/18/2008

09:48:47

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:49:32
Diagnostic Test	Pass	09:50:00
Air Blank	0.000	09:50:27
0.100 Cal Check	0.099	09:50:46
Air Blank	0.000	09:51:15
Subject Test	0.120*	09:54:24
Air Blank	0.000	09:54:55
Air Blank	0.000	09:55:23
Subject Test	XXX**	09:55:57
Air Blank	0.000	09:56:29

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

12/18/2008

09:56:50

Standard Lot#: 811201I

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:57:17
Cal Check	0.099	09:57:36
Air Blank	0.000	09:58:05

b of g

Cu

TS4

Clare Conley

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

758
Bh
12-24-08
Cu

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 12-18-08 TIME 0938
Intoxilyzer Serial # 80-000334 Location DUI HQ ODW

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
12/18/2008 09:38:12

Standard Lot#: 811201I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:38:39
Cal Check	0.099	09:38:58
Air Blank	0.000	09:39:28

8578
BY CU

Clare Conley
Operator's Signature