



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
8.31.2009

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 08/20/09 Cal 1 0.099 AC Cal 2 0.097 AC

Between the dates of 07/23/09 and 08/20/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


08/31/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 31, day of August, 2009.

A.M. Hernandez

Notary Public

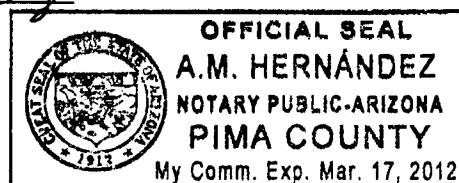

\*Includes documents specified in A.A.C. R13-10-104.A.5

187  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

134  
8.31.09

QA SPECIALIST Cheryl Straub AGENCY Tucson Police Department

DATE 08/20/09 TIME 12:18

INTOXILYZER SERIAL # 80-000334 LOCATION DUI Squad HQ

( ☒ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ☒ ) 1. Clock time check.
- ( ☒ ) 2. Date check.

OPERATIONAL TESTS

- ( ☒ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ☒ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ☒ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

\_\_\_\_\_

\_\_\_\_\_

\_\_\_\_\_

\_\_\_\_\_

SIGNATURE Cheryl Straub

287  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

08/20/2009

12:36:34

BM

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:37:41
Diagnostic Test	Pass	12:38:08
Air Blank	0.000	12:38:35
0.100 Cal Check	0.098	12:38:54
Air Blank	0.000	12:39:24
Subject Test	0.000	12:39:47
Air Blank	ABT*	12:39:56
Air Blank	0.000	12:40:24

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department*

387  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

08/20/2009

12:40:44

BY

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	12:41:33
Air Blank	IPS**	12:41:42

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department*

497  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

08/20/2009

12:42:05

154

Standard Lot#: 826001I

Last Changed By: C. CONLEY

#34517

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:42:58
Diagnostic Test	Pass	12:43:26
Air Blank	0.000	12:43:53
0.100 Cal Check	0.098	12:44:12
Air Blank	0.000	12:44:41
Subject Test	0.185*	12:47:50
Air Blank	0.000	12:48:22
Air Blank	0.000	12:48:50
Subject Test	XXX**	12:49:04
Air Blank	0.000	12:49:38

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department*

577cs

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

08/20/2009

12:50:06

AM

-----  
Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:50:34
Cal Check	0.097	12:50:52
Air Blank	0.000	12:51:22



\_\_\_\_\_  
Operator's Signature

647  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-16-104(A)

## ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

## STANDARD CALIBRATION CHECK PROCEDURE

B6  
8.31.09QA Specialist Cheryl STRAUB AGENCY Tucson Police DepartmentDATE 08/20/09 TIME 12:18Intoxilyzer Serial # 80-000334 Location Dut Squad HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

787  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

08/20/2009

12:18:48

154


-----  
Standard Lot#: 826001I

Last Changed By: C. CONLEY #34517  
-----

QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:19:16
Cal Check	0.099	12:19:34
Air Blank	0.000	12:20:04

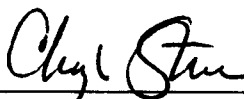

\_\_\_\_\_  
Operator's Signature