



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

By  
12-17-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/10/09 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 11/12/09 and 12/10/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

12/17/09

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17<sup>th</sup>, day of December, 2009.

A.M. Hernandez

Notary Public



\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY  
12-17-09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 12/10/09 TIME 1620

INTOXILYZER SERIAL # 80-000334 LOCATION Vehicle # 868-C15

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.

(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 8sc

184

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-100334  
12/10/2005 16:22:11

Standard Lot#: 842845  
Last Changed By: S. CIONE #50616

OPS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/21/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:23:02
Diagnostic Test	Pass	16:23:31
Air Blank	0.000	16:23:58
0.100 Cal Check	0.101	16:24:17
Air Blank	0.000	16:24:46
Subject Test	0.000	16:25:01
Air Blank	0.000	16:25:28
Five Minute Wait	AST*	
Air Blank	0.000	16:26:09

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 8 sc

BA

INTOXILYZER 8500  
Location: TUCSON POLICE DEPT  
Serial Number: 60-100334  
12/15/2009 16:26:15

Standard Lot#: 642645  
Last Changed By: S. CIONE #50616

Off: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SOAP\_SOP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RF1*	16:26:44
Air Blank	1PS**	16:26:52

\*RF1 Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 8 SC

TS4

INITIALIZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 60-000334  
12/10/2009 16:26:58

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

ORG: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SOAP,SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Depriuation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:27:44
Diagnostic Test	Fail*	16:28:14
Air Blank	0.000	16:28:43

\* Temperature Regulation Test

Not a Successfully  
Completed  
Test Sequence

-windows in vehicle  
were down + windy  
outside. Rolled up  
windows after this  
Diagnostic Test failure.

SC

5 of 8 sc

AS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
12/10/2009 16:28:57

Standard Lock: 642845  
Last Changed By: S. CIONE #50616

ORS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SORP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Denaturation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:29:44
Diagnostic Test	Pass	16:30:12
Air Blank	0.000	16:30:40
0.100 Cal Check	0.100	16:30:50
Air Blank	0.000	16:31:20
Subject Test	0.123*	16:34:36
Air Blank	0.000	16:35:09
Air Blank	0.000	16:35:36
Subject Test	XXX**	16:35:00
Air Blank	0.000	16:36:32

\*Deficient Sample  
\*\*Invalid Sample  
Begin new denaturation period

Not a Successfully  
Completed  
Test Sequence

6 of 8 sc

TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 00-000334  
12/15/2009 16:37:34

Standard Lot#: 642645  
Last Changed By: S. CIONE #50616

OPR: S. CIONE #50616  
TUCSON PD

Test	gr210L	Time
Air Blank	0.000	16:38:01
Cal Check	0.099	16:38:20
Air Blank	0.000	16:38:50

*Sarah Cione*  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

7 of 8<sub>5c</sub>

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

By  
12-17-09

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 12/10/09 TIME 1620  
Intoxilyzer Serial # 80-000334 Location Vehicle # 868 - C15

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

Page 8 of 50

184

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
12/10/2009 16:20:49

Standard Lot#: 642645  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:21:16
Cal Check	0.099	16:21:35
Air Blank	0.000	16:22:04

*Sarah Cione*  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

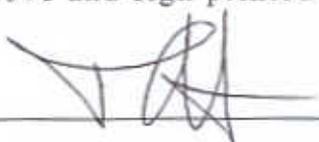
~~1016?~~  
1 of 2 sc

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

TSY  
12-17-09

Date: 11-17-09  
 QAS: O'HARA # 38876  
 Intoxilyzer Serial Number: 80-000334  
 New Lot Number: 826001I  
 New Expiration Date: 25 SEP 2010

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
 Select "D" for dry gas and enter the following information:  
     Enter Target Value: 0.100  
     Enter Control Lot #: From the tank label  
     Enter Exp Date: From the tank label (mm/dd/yyyy)  
     Edit info Y/N? Y or N  
     Update Std Info? Y ENTER  
     "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
 Select "D" for dry control test. After selection is made press  
 ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 101 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE 

2 of 25c

154

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
11/17/2009 23:16:40

Standard Lot#: 826011  
Last Changed By: T. O HARR #38876

OPS: T. O HARR #38876  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	23:17:09
Cal Check	0.101	23:17:28
Air Blank	0.000	23:17:58

  
Operator's Signature 38876



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 2 sc

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

TS4  
12-17-09

Date: 11/20/09

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000334

New Lot Number: 642845

New Expiration Date: 7/23/2011

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:  
 Enter Target Value: 0.100  
 Enter Control Lot #: From the tank label  
 Enter Exp Date: From the tank label (mm/dd/yyyy)  
 Edit info Y/N? Y or N  
 Update Std Info? Y ENTER  
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

2 of 2 sc

154

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 92-200334  
11/20/2005 10:53:29

Standard Lots: 642945  
Last Changed By: S. CIONE #50616

QRS: S. CIONE #50616  
TUCSON PD

Test	g/210L	Time
Air Blank	0.300	10:53:56
Cal Check	0.100	10:54:15
Air Blank	0.800	10:54:45

*Sarah Cione*  
Operator's Signature