



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

Tsy  
10-18-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 87-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10-14-10 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 9-16-10 and 10-14-10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


10-15-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 15, day of October, 2010.

Laura B Hodgson

Notary Public

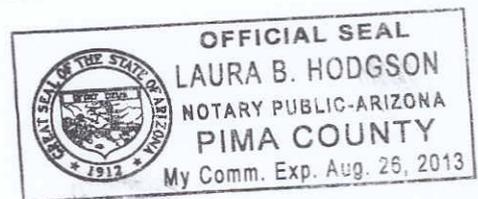

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1078  
lu  
TSK 10-18-2010

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 10-14-10 TIME 0830

INTOXILYZER SERIAL # 80-000334 LOCATION ODW DM HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
10/14/2010 08:32:35

258  
EW  
184

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:33:24
Diagnostic Test	Pass	08:33:51
Air Blank	0.000	08:34:19
0.100 Cal Check	0.099	08:34:38
Air Blank	0.000	08:35:07
Subject Test	0.000	08:35:40
Air Blank	ABT*	08:35:46
Air Blank	0.000	08:36:14

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
10/14/2010

08:37:23

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:38:07
Air Blank	IPS**	08:38:16

\*RFI Detect  
\*\*Improper Sample

358  
OK  
B4

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
10/14/2010 08:39:38

4678  
CW  
B4

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:40:38
Diagnostic Test	Pass	08:41:06
Air Blank	0.000	08:41:33
0.100 Cal Check	0.098	08:41:52
Air Blank	0.000	08:42:21
Subject Test	0.107*	08:45:30
Air Blank	IPS**	08:45:37
Air Blank	0.000	08:46:08

\*Deficient Sample  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000334  
10/14/2010 08:46:31

*508  
lu*

Standard Lot#: 642845  
Last Changed By: S. CIONE #50616

*Bh*

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:47:25
Diagnostic Test	Pass	08:47:52
Air Blank	0.000	08:48:20
0.100 Cal Check	0.098	08:48:38
Air Blank	0.000	08:49:08
Subject Test	XXX*	08:49:38
Air Blank	0.000	08:50:11

\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
10/14/2010 09:11:58

bbs  
cu  
134

-----  
Standard Lot#: 642845  
Last Changed By: S. CIONE #50616  
-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:12:26
Cal Check	0.098	09:12:44
Air Blank	0.000	09:13:14

Clare Conley  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

758  
Cu

TSY  
10-18-2010

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 10-14-10 TIME 0830  
Intoxilyzer Serial # 80-000334 Location ODW DWI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000334  
10/14/2010 08:30:41

8/18  
CU

-----  
Standard Lot#: 642845  
Last Changed By: S. CIONE #50616  
-----

154

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:31:09
Cal Check	0.099	08:31:27
Air Blank	0.000	08:31:57

*Clare Conley*

\_\_\_\_\_  
Operator's Signature