



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
8.22.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-18-11 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 7-21-11 and 8-18-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.

The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


8-18-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 18th, day of August, 20 11.

Susan A. Hitchye

Notary Public

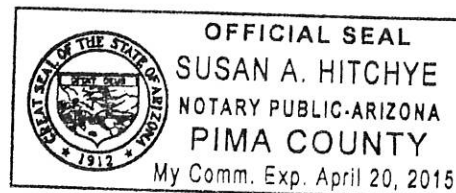

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

67
Cu
B4
8.22.2011

QA SPECIALIST Mare Conley AGENCY Tucson Police Department

DATE 8-18-11 TIME 0851

INTOXILYZER SERIAL # 80-000334 LOCATION ODW DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS None.

SIGNATURE Mare Conley

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

08/18/2011

08:55:44

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:56:37 |
| Diagnostic Test | Pass | 08:57:05 |
| Air Blank | 0.000 | 08:57:32 |
| 0.100 Cal Check | 0.098 | 08:57:51 |
| Air Blank | 0.000 | 08:58:20 |
| Subject Test | 0.000 | 08:58:52 |
| Air Blank | ABT* | 08:58:58 |
| Air Blank | 0.000 | 08:59:27 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
08/18/2011

08:59:36

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 09:00:23 |
| Diagnostic Test | Pass | 09:00:51 |
| Air Blank | RFI* | 09:00:58 |
| Air Blank | IPS** | 09:01:08 |

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-000334

08/18/2011

09:01:30

Standard Lot#: 20910100A1

Last Changed By: S. CIONE

#50616

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 09:02:20 |
| Diagnostic Test | Pass | 09:02:48 |
| Air Blank | 0.000 | 09:03:15 |
| 0.100 Cal Check | 0.098 | 09:03:34 |
| Air Blank | 0.000 | 09:04:04 |
| Subject Test | 0.053* | 09:07:12 |
| Air Blank | 0.000 | 09:07:43 |
| Air Blank | 0.000 | 09:08:11 |
| Subject Test | XXX** | 09:08:32 |
| Air Blank | 0.000 | 09:09:05 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

08/18/2011

09:09:20


Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 09:09:47 |
| Cal Check | 0.098 | 09:10:06 |
| Air Blank | 0.000 | 09:10:36 |

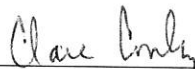

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 8-18-11 TIME 0851
Intoxilyzer Serial # 80-000334 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000334

08/18/2011

08:51:44


Standard Lot#: 20910100A1

Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 08:52:12 |
| Cal Check | 0.098 | 08:52:30 |
| Air Blank | 0.000 | 08:53:00 |

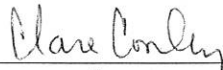

Operator's Signature