



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
3.20.2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/21/12 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 2/2/12 and 2/21/12 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on 2/21/12.
 The instrument was returned to service on _____.

Sarah Cione
 Sarah L. Cione

3/9/12
 ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this *9th* day of *March*, 20*12*.

Elma Escarcega
 Notary Public

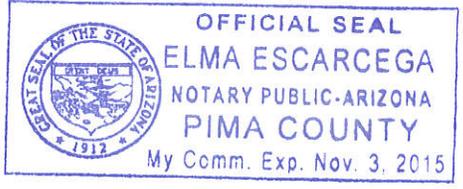

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 9_{5c}

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4
3.20.2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 2/21/12 TIME 1419
INTOXILYZER SERIAL # 80-000334 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 9 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
02/21/2012

14:21:42

BSH

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:22:24
Diagnostic Test	Pass	14:22:52
Air Blank	0.000	14:23:19
0.100 Cal Check	0.099	14:23:38
Air Blank	0.000	14:24:08
Subject Test	0.000	14:24:29
Air Blank	0.000	14:24:56
Five Minute Wait	ABT*	
Air Blank	0.000	14:25:30

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 9 sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
02/21/2012 14:25:35

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:26:10
Air Blank	IPS**	14:26:17

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 95c

BSH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000334
02/21/2012 14:26:30

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:27:15
Diagnostic Test	Pass	14:27:43
Air Blank	0.000	14:28:10
0.100 Cal Check	0.099	14:28:29
Air Blank	0.000	14:28:58
Subject Test	0.093*	14:32:07
Air Blank	0.000	14:32:39
Air Blank	0.000	14:33:06
Subject Test	XXX**	14:33:23
Air Blank	0.000	14:33:58

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 9 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
02/21/2012 14:34:28

TS4

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:34:55
Cal Check	0.099	14:35:14
Air Blank	0.000	14:35:43

Sarah Cione

Operator's Signature

EXHIBIT G-2

6 of 9 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
3.20.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 2/21/12 TIME 1419
Intoxilyzer Serial # 80-000334 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
02/21/2012 14:19:41

7 of 9 SC
Tsh

Standard Lot#: 20910100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:20:08
Cal Check	0.100	14:20:27
Air Blank	0.000	14:20:57

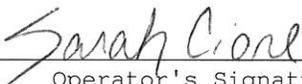

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9 sc

INTOXILYZER SERVICE RECORD

Bh
3.20.2012

Date: 2/21/12

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-000334

Reason For Service: Intox had a couple of modem failures. Notified by officer Marin on 2/17/12.

Date Out of Service: 2/21/12

Calibration Check completed prior to removal from service

Date: 2/21/12

SQAP completed prior to removal from service

Date: 2/21/12

Corrective Action: Sent to CMI, Inc.



Service Evaluation Form

9 of 9 sc

TS4
3.20.2012

This form **MUST** be completed and enclosed with instrument to be serviced.
Failure to complete and return this form may cause delays in service.

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)
 Name Terry Gallegos Phone: (520) 837-7677
 Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:
Finance / Accounts Payable Tucson Police Department Crime Lab
Po Box 27450 1306 W. Miracle Mile
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 80-000334 Instrument Model: I8000

4. Detailed Description of Problem:
Periodic modem failures during Diagnostics. Also updated software to 8105.48,
but then wouldn't pass Analytical Test.

*****Hazardous Material Warning! - DO NOT return gas cylinders with instrument!*****

I Authorize Repairs Up To: All \$250 \$500 \$750 Other \$ _____
 Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:
Terry Gallegos Toxicology Coordinator
 Name (Please Print) Title
Terry Gallegos 3.1.2012
 Signature Date

Ship to:
 CMI, Inc.
 Attn: Service Dept.
 316 East Ninth Street
 Owensboro, KY 42303

No, please send estimate before repairs are made.
 Note: An estimate will be faxed before performing any repairs and may cause delays in service.
 An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.