



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
8-27-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000334 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/16/12 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 7/19/12 and 8/16/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


8/21/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 21st, day of August, 2012.

Edna Escarcega

Notary Public

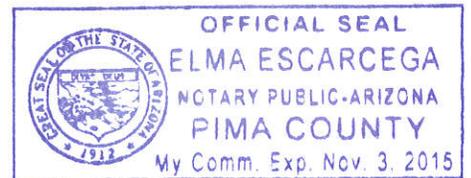

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8 sc
TSg
8-27-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 8/16/12 TIME 1102

INTOXILYZER SERIAL # 80-000334 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 8_{sc}
Bh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
08/16/2012 11:05:43

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:06:26
Diagnostic Test	Pass	11:07:00
Air Blank	0.000	11:07:28
0.100 Cal Check	0.100	11:07:48
Air Blank	0.000	11:08:18
Subject Test	0.000	11:08:46
Air Blank	0.000	11:09:14
Wait	ABT*	
Air Blank	0.000	11:09:49

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8
SC
Bh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
08/16/2012 11:18:30

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:18:58
Air Blank	IPS**	11:19:06

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 sc
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
08/16/2012 11:21:01

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:21:46
Diagnostic Test	Pass	11:22:20
Air Blank	0.000	11:22:47
0.100 Cal Check	0.101	11:23:08
Air Blank	0.000	11:23:38
Subject Test	XXX*	11:24:19
Air Blank	0.000	11:24:53

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 8 sc
TJ4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
08/16/2012 11:28:05

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:28:48
Diagnostic Test	Pass	11:29:22
Air Blank	0.000	11:29:50
0.100 Cal Check	0.100	11:30:11
Air Blank	0.000	11:30:40
Subject Test	0.046*	11:33:48
Air Blank	0.000	11:34:19
Air Blank	0.000	11:34:46
Subject Test	ABT**	11:35:17
Air Blank	0.000	11:35:45

*Deficient Sample
**Sequence Aborted

Not a Successfully Completed Test Sequence

6 of 85c
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
08/16/2012 11:36:40

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:37:09
0.100 Cal Check	0.100	11:37:29
Air Blank	0.000	11:37:59

Preventative Maintenance
Performed -- Timer Reset

7 of 8 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

TS
8-27-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 8/16/12 TIME 1102
Intoxilyzer Serial # 80-000334 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
08/16/2012 11:02:32

8 of 8 sc
TH

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:03:01
0.100 Cal Check	0.100	11:03:21
Air Blank	0.000	11:03:51

Preventative Maintenance
Performed -- Timer Reset