



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

By
10-29-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-0003734 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/11/12 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 9/13/12 and 10/11/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione



10/18/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 18th day of October, 20 12.

Elma Escarcega

Notary Public

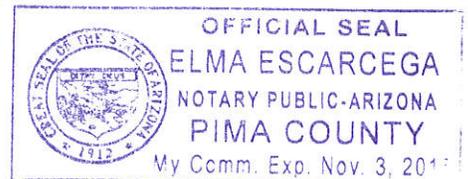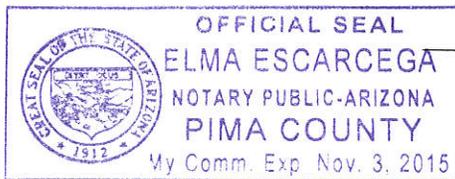

*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7_{sc}

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TS4
10-29-2012

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 10/11/12 TIME 1058

INTOXILYZER SERIAL # 80-000334 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 75c
Bh

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
10/11/2012 11:02:33

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:03:15
Diagnostic Test	Pass	11:03:49
Air Blank	0.000	11:04:16
0.100 Cal Check	0.100	11:04:37
Air Blank	0.000	11:05:06
Subject Test	0.000	11:05:19
Air Blank	0.000	11:05:47
Wait	ABT*	
Air Blank	0.000	11:06:51

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}

Sh

INTOXILYZER 8000
 Location: TUCSON POLICE DEPT
 Serial Number: 80-000334
 Core Version: 8105.48
 10/11/2012 11:06:58

Standard Lot#: 23511100A2
 Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
 TUCSON PD

Operator: S. CIONE #50616
 TUCSON PD

Subject: SQAP, SQAP
 DOB: 11/11/1980
 Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:07:46
Diagnostic Test	Pass	11:08:20
Air Blank	0.000	11:08:47
0.100 Cal Check	RFI*	11:09:06
Air Blank	IPS**	11:09:13

*RFI Detect
 **Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc
TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
10/11/2012 11:17:37

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:18:27
Diagnostic Test	Pass	11:19:01
Air Blank	0.000	11:19:28
0.100 Cal Check	0.100	11:19:49
Air Blank	0.000	11:20:18
Subject Test	0.087*	11:23:26
Air Blank	0.000	11:23:58
Air Blank	0.000	11:24:25
Subject Test	XXX**	11:24:47
Air Blank	0.000	11:25:23

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
10/11/2012 11:28:13

5 of 7 sc
Bh

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:28:41
0.100 Cal Check	0.100	11:29:02
Air Blank	0.000	11:29:31

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 50
TS4
10-29-2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 10/11/12 TIME 1058
Intoxilyzer Serial # 80-000334 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 75c

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-000334
Core Version: 8105.48
10/11/2012 10:58:31

Standard Lot#: 23511100A2
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:00:30
0.100 Cal Check	0.101	11:00:51
Air Blank	0.000	11:01:20

Preventative Maintenance
Performed -- Timer Reset