



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program. B4

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.


Date: 12/1/06 Cal 1 0.099 AC Cal 2 0.099 AC

The instrument was placed into service on 12/1/06.

Jennifer L. Coombs  
Jennifer L. Coombs  
12/27/06

ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 27, day of December, 2006.

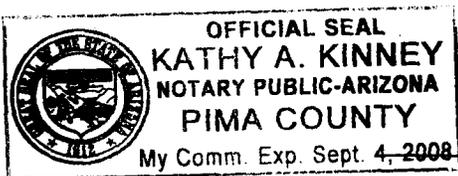

Kathy A. Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 8  
JL

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE <sup>DU</sup>

QA SPECIALIST Jennifer Coombs 46160 AGENCY Tucson Police Department  
DATE December 1, 2006 TIME 1312  
INTOXILYZER SERIAL # 80-000712 LOCATION LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Jennifer L. Coombs

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
12/01/2006

13:15:47

ph

2 of 8  
file

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:16:45
Diagnostic Test	Pass	13:17:12
Air Blank	0.000	13:17:39
0.100 Cal Check	0.099	13:17:58
Air Blank	0.000	13:18:28
Subject Test	0.000	13:18:45
Air Blank	0.000	13:19:12
Five Minute Wait	ABT*	
Air Blank	0.000	13:19:59

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
12/01/2006

13:20:55

DH

3 of 8  
file

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:21:47
Diagnostic Test	Pass	13:22:14
Air Blank	IPS*	13:22:22
Air Blank	0.000	13:22:50

\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
12/01/2006

13:22:56

bh

428  
file

---

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

---

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:23:47
Diagnostic Test	Pass	13:24:14
Air Blank	0.000	13:24:41
0.100 Cal Check	0.100	13:25:00
Air Blank	0.000	13:25:30
Subject Test	DEF*	13:28:38
Air Blank	0.000	13:29:06
Air Blank	0.000	13:29:33
Subject Test	RFI**	13:30:00
Air Blank	0.000	13:30:28

\*Deficient Sample

\*\*RFI Detect

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-000712  
12/01/2006

13:32:39

PH

5 of 8  
file

Standard Lot#: 502602I  
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160  
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160  
TUCSON PD CRIME LAB

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:34:01
Diagnostic Test	Pass	13:34:29
Air Blank	0.000	13:34:56
0.100 Cal Check	0.099	13:35:15
Air Blank	0.000	13:35:44
Subject Test	XXX*	13:36:12
Air Blank	0.000	13:36:45

\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

12/01/2006

13:37:07

6 of 8  
file  
m

-----  
Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160  
-----

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	13:37:34
Cal Check	0.099	13:37:53
Air Blank	0.000	13:38:22

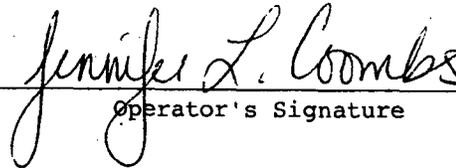
  
\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

7 of 8  
me

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

ph

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Jennifer Coombs 46660 AGENCY Tucson Police Department  
DATE December 1, 2006 TIME 1312  
Intoxilyzer Serial # 80-000712 Location LAB

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Jennifer L. Coombs

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

12/01/2006

13:12:02

124

8 of 8  
JL

-----  
Standard Lot#: 502602I

Last Changed By: J. COOMBS #46160

-----  
QAS: J. COOMBS #46160

TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	13:12:30
Cal Check	0.099	13:12:49
Air Blank	0.000	13:13:18

  
\_\_\_\_\_  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

1 of 1  
DH  
JL

Intoxilyzer Serial Number: 80-000712

Initial Service Date: 12-1-06

CMI, Inc. Certificate of Calibration Date: 4-14-04