



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

WMS
5-25-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4.26.2007 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 3.26.2007 and 4.26.2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S. Gallegos
Terry S. Gallegos

5.24.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 24, day of May, 2007.

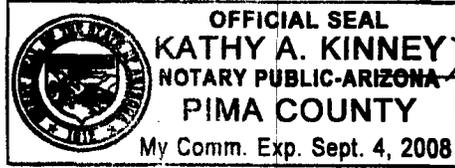

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

1 of 8
TS4
5-25-07

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 4-26-2007 TIME 1557

INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
04/26/2007

16:00:20

2 of 8
B4 YMB
5-2507

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:01:08
Diagnostic Test	Pass	16:01:36
Air Blank	0.000	16:02:03
0.100 Cal Check	0.100	16:02:22
Air Blank	0.000	16:02:51
Subject Test	0.000	16:03:34
Air Blank	0.000	16:04:01
Five Minute Wait	ABT*	
Air Blank	0.000	16:04:41

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
04/26/2007

16:05:17

30 of 8
5-2507
TS4

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:06:06
Diagnostic Test	Pass	16:06:34
Air Blank	0.000	16:07:01
0.100 Cal Check	RFI*	16:07:13
Air Blank	IPS**	16:07:28

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
04/26/2007

16:07:42

4 of 8
754
32507

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:08:31
Diagnostic Test	Pass	16:08:58
Air Blank	0.000	16:09:26
0.100 Cal Check	0.100	16:09:44
Air Blank	0.000	16:10:14
Subject Test	0.137*	16:13:22
Air Blank	0.000	16:13:54
Air Blank	0.000	16:14:21
Subject Test	0.414	16:16:59
Air Blank	0.000	16:17:37
Five Minute Wait	ABT**	
Air Blank	0.000	16:18:35

*Deficient Sample
**Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-000712
04/26/2007

16:19:30

5 of 8
TS4
5/26/07

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:20:21
Diagnostic Test	Pass	16:20:48
Air Blank	0.000	16:21:16
0.100 Cal Check	0.101	16:21:35
Air Blank	0.000	16:22:04
Subject Test	XXX*	16:22:24
Air Blank	0.000	16:22:58

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 8
TS4 5-25-07

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 4.26.2007 TIME 1557
Intoxilyzer Serial # 80-000712 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-000712

04/26/2007

15:57:41

8078
TS/ 52507

Standard Lot#: 635301I

Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test g/210L Time

Air Blank 0.000 15:58:08
Cal Check 0.100 15:58:27
Air Blank 0.000 15:58:56



Operator's Signature