



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

CU  
10.26.07

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9.19.2007 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 9.5.2007 and 9.19.2007 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.


Terry S Gallegos

Terry S. Gallegos

10.25.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 25, day of October, 2007.

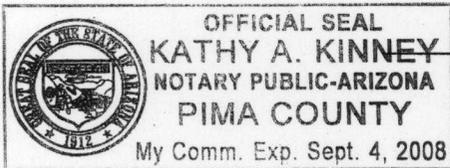

Kathy A Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

*CU*  
*10.26.07*

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 9.19.2007 TIME 1756

INTOXILYZER SERIAL # 80-000712 LOCATION DUI HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓✓
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE *Terry Gallegos*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
09/19/2007 18:03:53

-----  
Standard Lot#: 6353011  
Last Changed By: G. EPPLEY #26227  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:05:14
Diagnostic Test	Pass	18:05:43
Air Blank	0.000	18:06:11
0.100 Cal Check	0.098	18:06:29
Air Blank	0.000	18:06:59
Subject Test	0.000	18:08:09
Air Blank	0.000	18:08:36
Five Minute Wait	ABT*	
Air Blank	0.000	18:09:11

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

CU

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
09/19/2007 18:09:31

-----  
Standard Lot#: 6353011  
Last Changed By: G. EPPLEY #26227  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	18:10:21
Diagnostic Test	Pass	18:10:50
Air Blank	0.000	18:11:18
0.100 Cal Check	0.098	18:11:36
Air Blank	0.000	18:12:06
Subject Test	RFI*	18:12:51
Air Blank	IPS**	18:13:11

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

CM

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
09/19/2007 18:14:11

-----  
Standard Lot#: 6353011  
Last Changed By: G. EPPLEY #26227  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	18:15:10
Diagnostic Test	Pass	18:15:38
Air Blank	0.000	18:16:06
0.100 Cal Check	0.099	18:16:24
Air Blank	0.000	18:16:54
Subject Test	0.137*	18:20:02
Air Blank	0.000	18:20:34
Air Blank	0.000	18:21:01
Subject Test	XXX**	18:21:18
Air Blank	0.000	18:21:52

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

CM

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
09/19/2007 18:23:56

-----  
Standard Lot#: 6353011  
Last Changed By: G. EPPLEY #26227  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:24:23
Cal Check	0.098	18:24:42
Air Blank	0.000	18:25:12

  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104 (A)

CM  
10-26-07


ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 20047 AGENCY Tucson Police Department  
DATE 9-19-2007 TIME 1756  
Intoxilyzer Serial # 80-000712 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 098 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

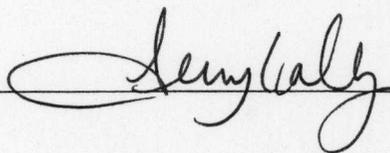SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

24

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-000712  
09/19/2007 17:56:12

-----  
Standard Lot#: 6353011  
Last Changed By: G. EPPLEY #26227  
-----

OAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	17:56:40
Cal Check	0.098	17:56:58
Air Blank	0.000	17:57:28

  
-----  
Operator's Signature